

TPC 3D




Software Installation

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

Download IntraoralScan_Aidite_setup.exe


Download Link : <https://service.aiditecloud.com/main>

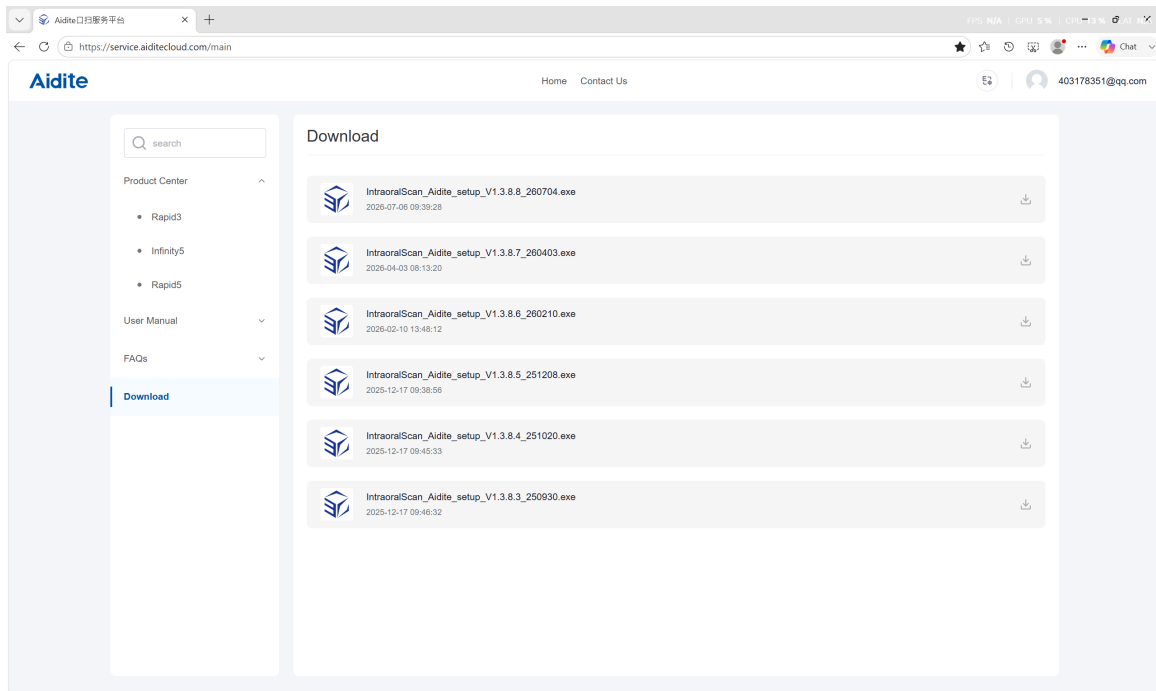


Alternatively, you may locate this file on the USB flash drive included in the original package.

Run IntraoralScan_Aidite_setup.exe

Select Language

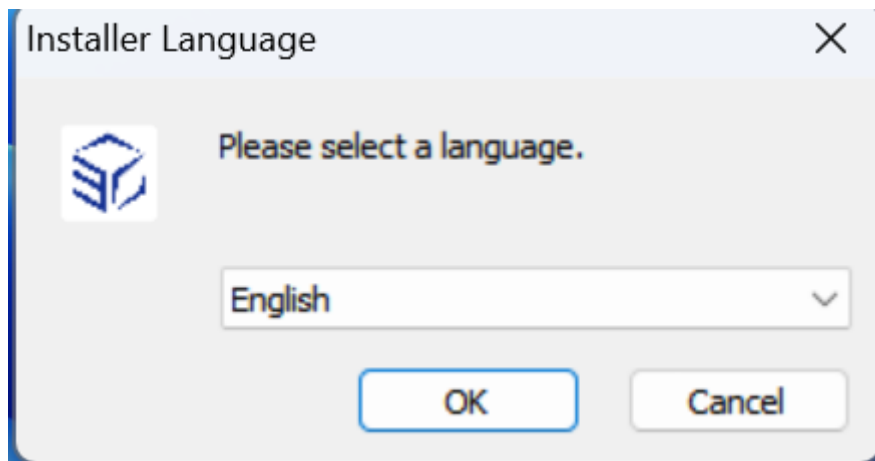


Agree to EULA

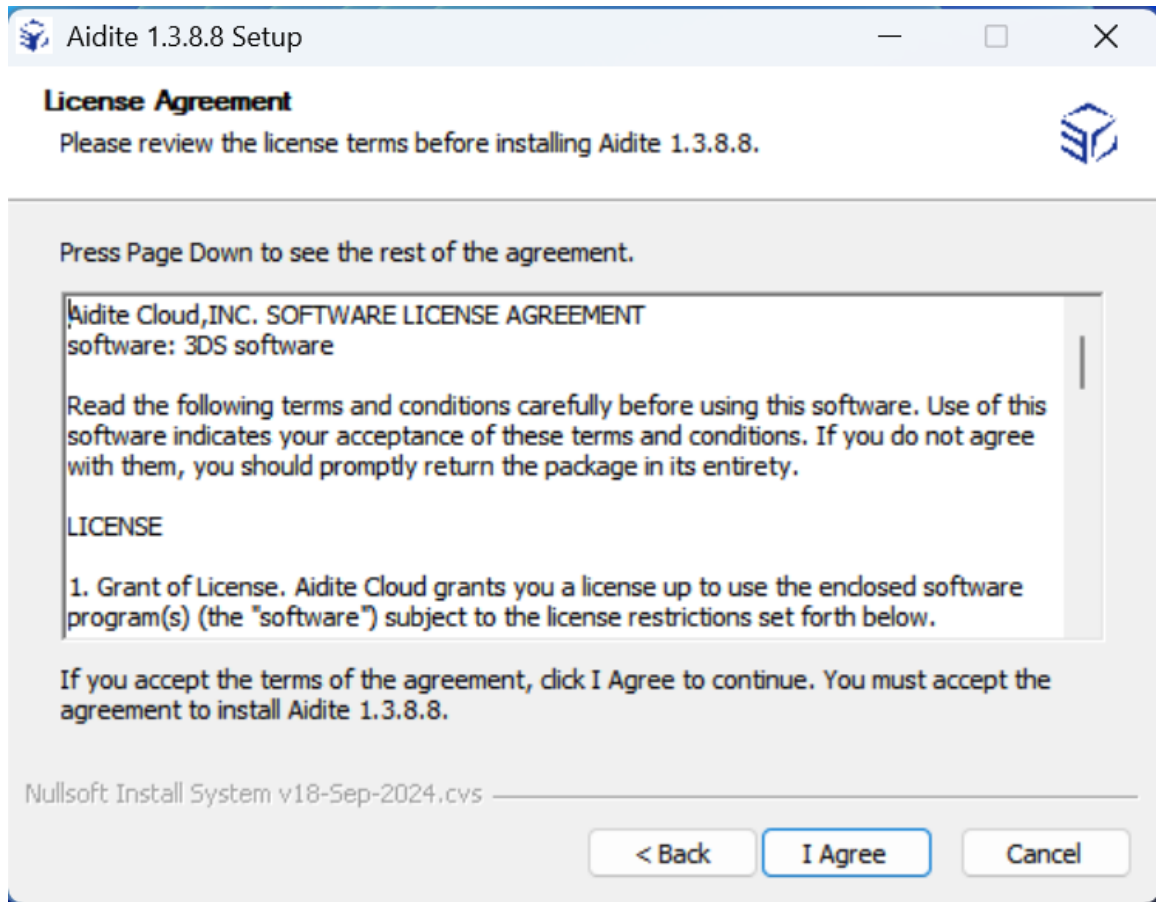


Start Installation

The program will install automatically and the process may take several minutes. Do not shut down your computer or unplug the power supply before installation completes.

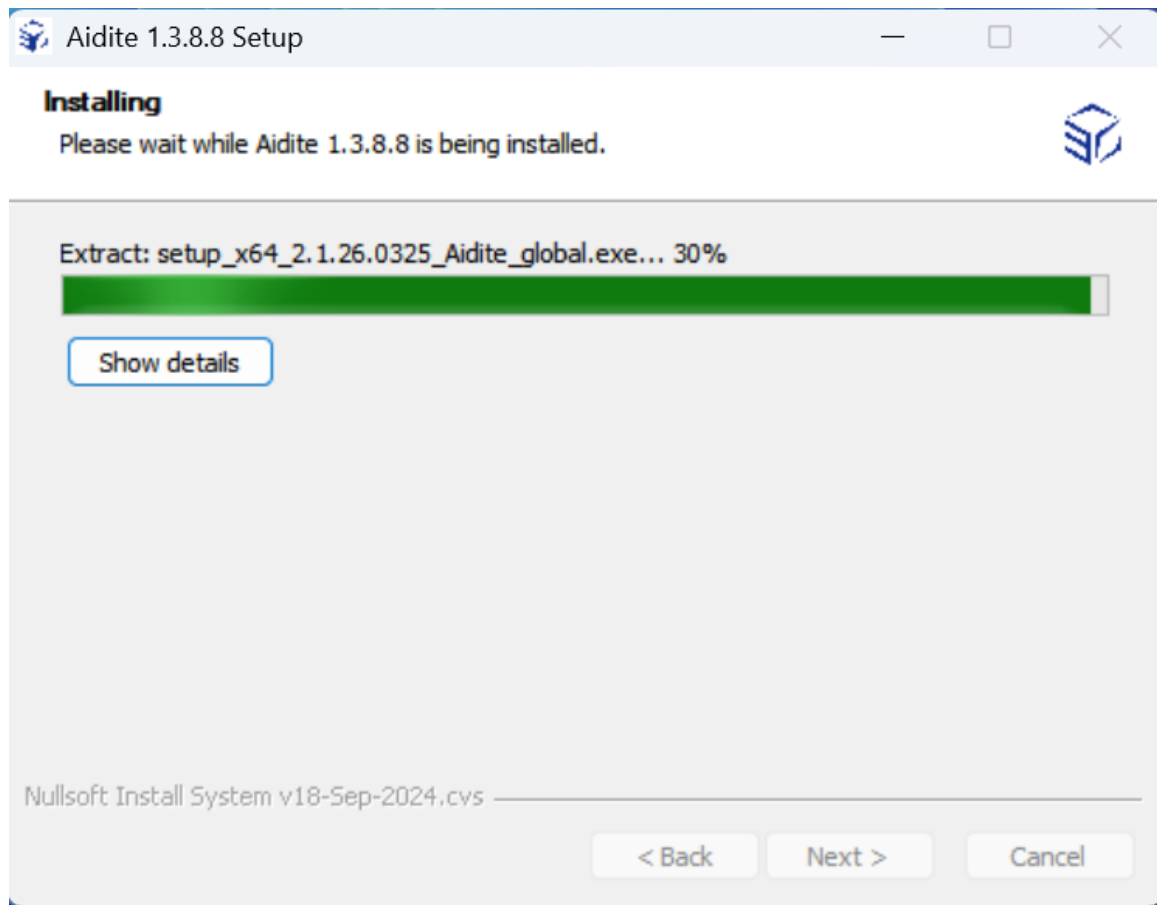



Installation Complete

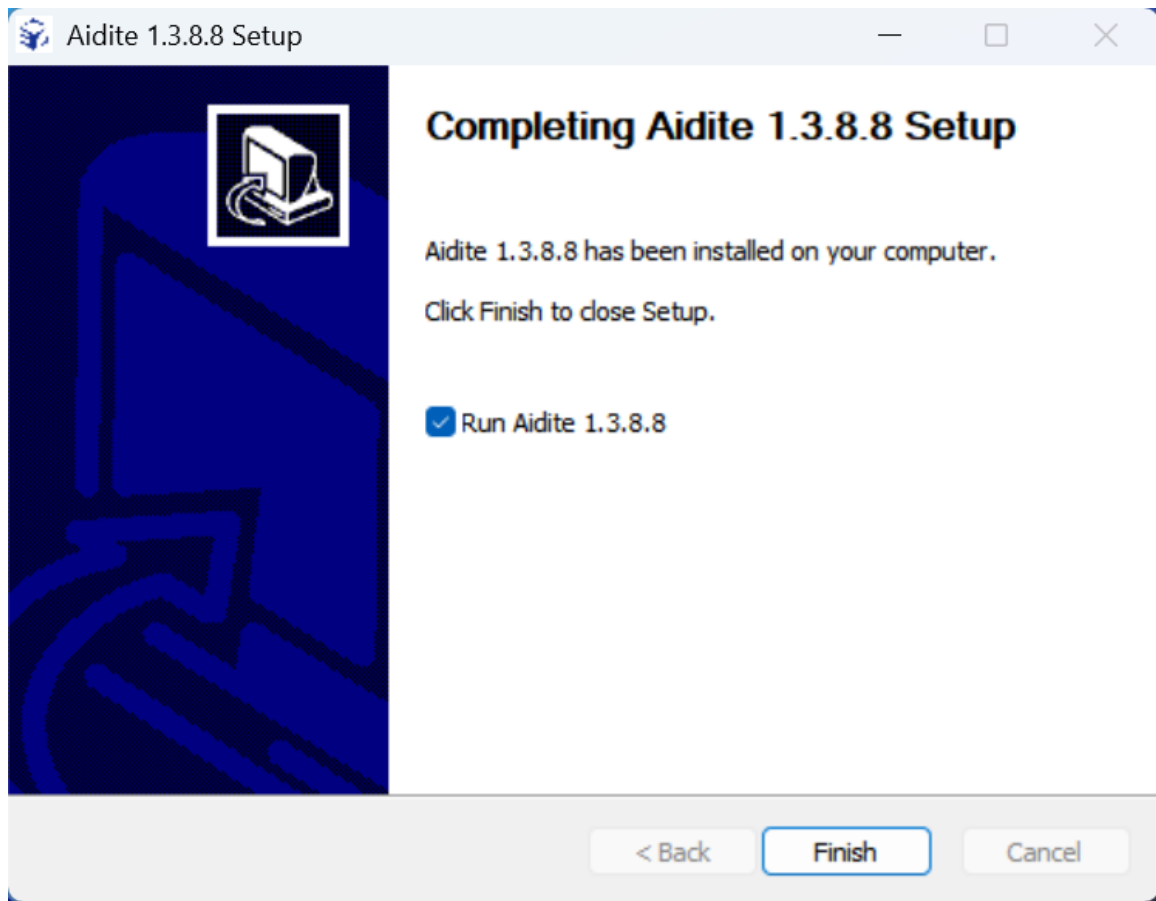

Software Update

You can check your software version and perform updates in the Settings menu on the software startup screen.

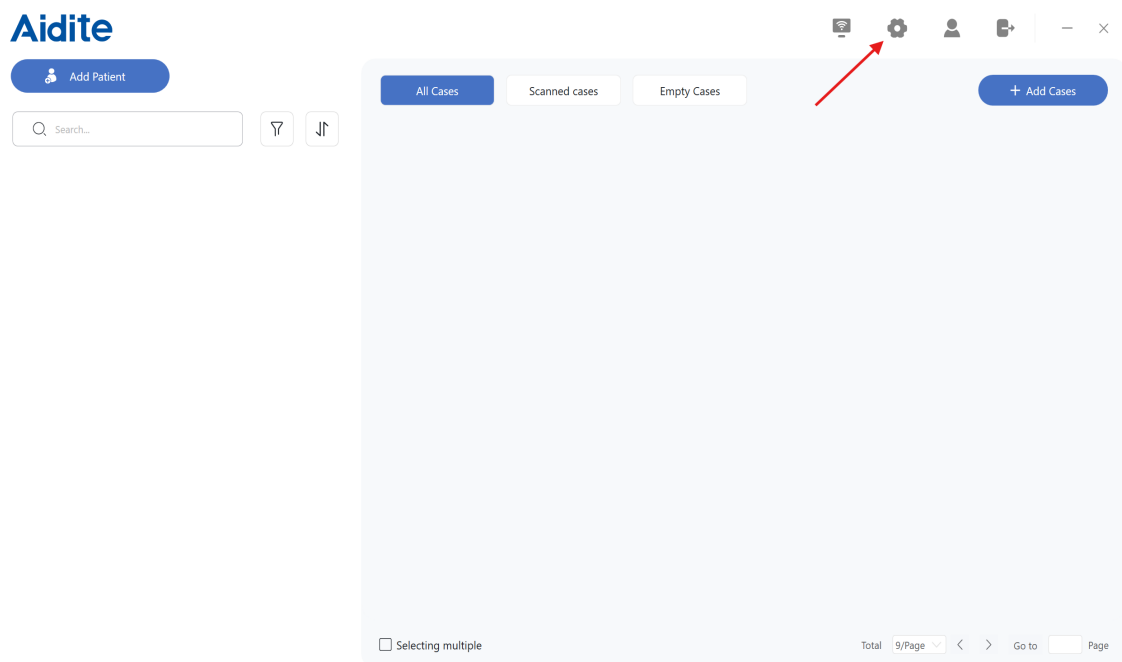

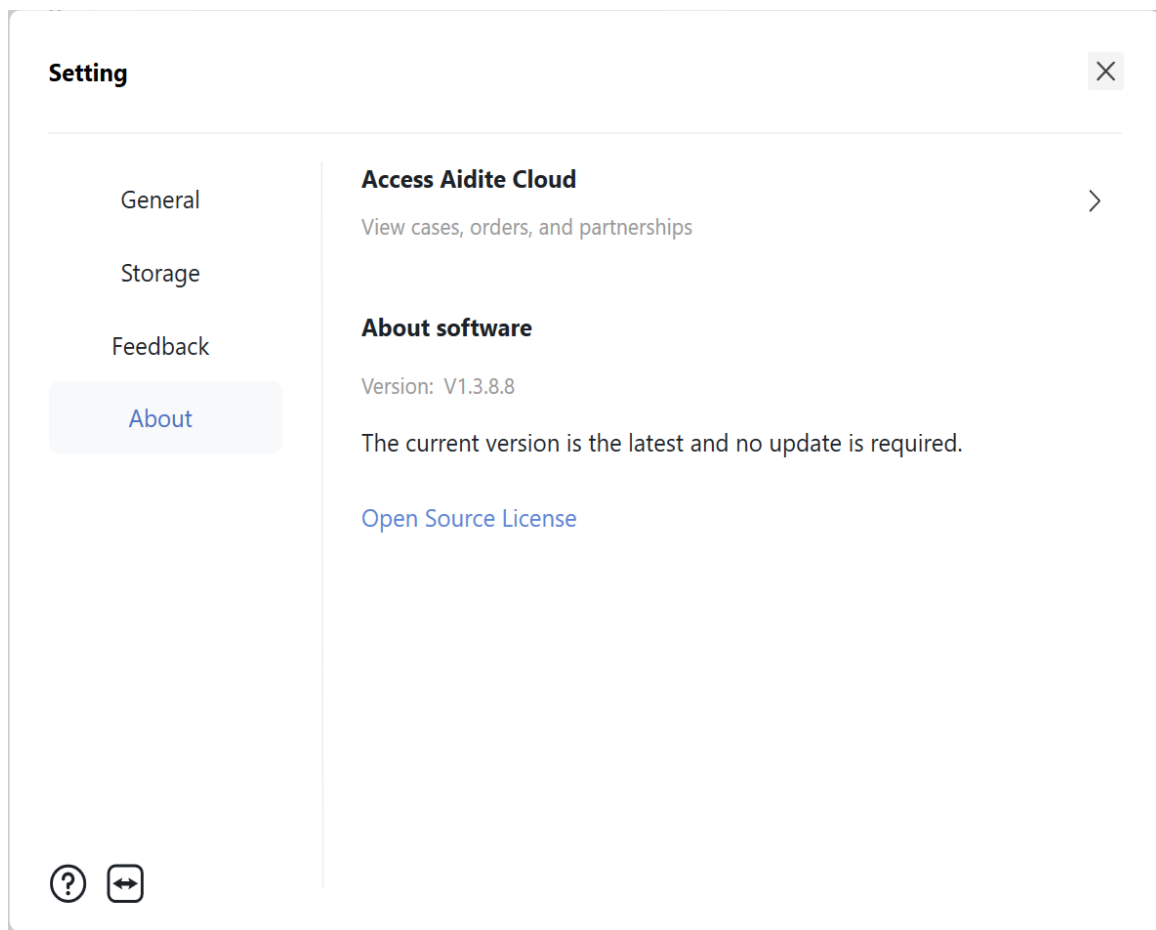

Alternatively, you may update the software by downloading the latest installation package from the website <https://service.aiditecloud.com/main> and reinstalling it. There is no need to worry about storage space, as our installer will automatically clean up the old software version.