

TPC 3D



Account Registration

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

Open the software and enter the login interface.



Click "Register" and enter the registration interface.

Aidite Oral Scanning

Welcome to Aidite Oral Scanning

☐ Remember Password

log in

Forgot Password **Register**

* By clicking on login, you agree to the terms of service and privacy policy

Or you can open <https://intraoralscan.aiditecloud.com/>

Dental Digital Center

Welcome to AICLOUD Platform

Log in

Forgot password **Register** 中/en

*By log in I consent to the [cloud data processing agreement](#)
Copyright © Aidite (Qinhuangdao) Technology Co.,Ltd | ICP Bei 17008635 No.5 V1.2.1

Select a role type: outpatient or factory, and then click "Next". Outpatients can send orders to factories which can also receive and download orders.

REGISTER

① Select character

Please select your character


Lab


Clinic

← Back to log in

Next step →

Select or fill in your country, organization name, and username. Click "Next".

REGISTER

② Fill in the basic information

Country

Please select your country

Clinic Name

Please enter your clinic name

User Name

Please enter your name

← Previous step

Next step →

Fill in your email and set a password. Click "Next".

REGISTER

③ Determine email password

E-mail

Please enter your email address

Password

Enter a password of at least 6 characters

Confirm Password

Please enter your password again

← Previous step

Next step →

The software will send a verification email to your registered email address.

REGISTER

Congratulations!

A confirmation email has been sent to:


Please check your email and follow the instructions to complete your registration and begin to use

Didn't receive a confirmation email?

[Click to resend](#)

← Back to log in

Open the email and verify your identity to complete the account registration.

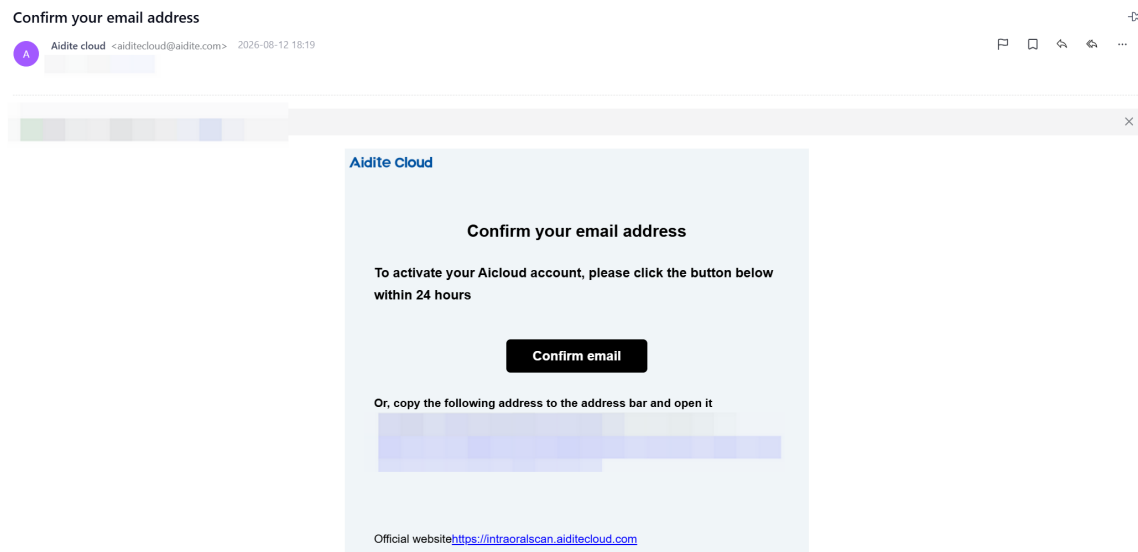

Return to the login interface and enter your account and password.