


TPC 3D



Scanner Installation

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

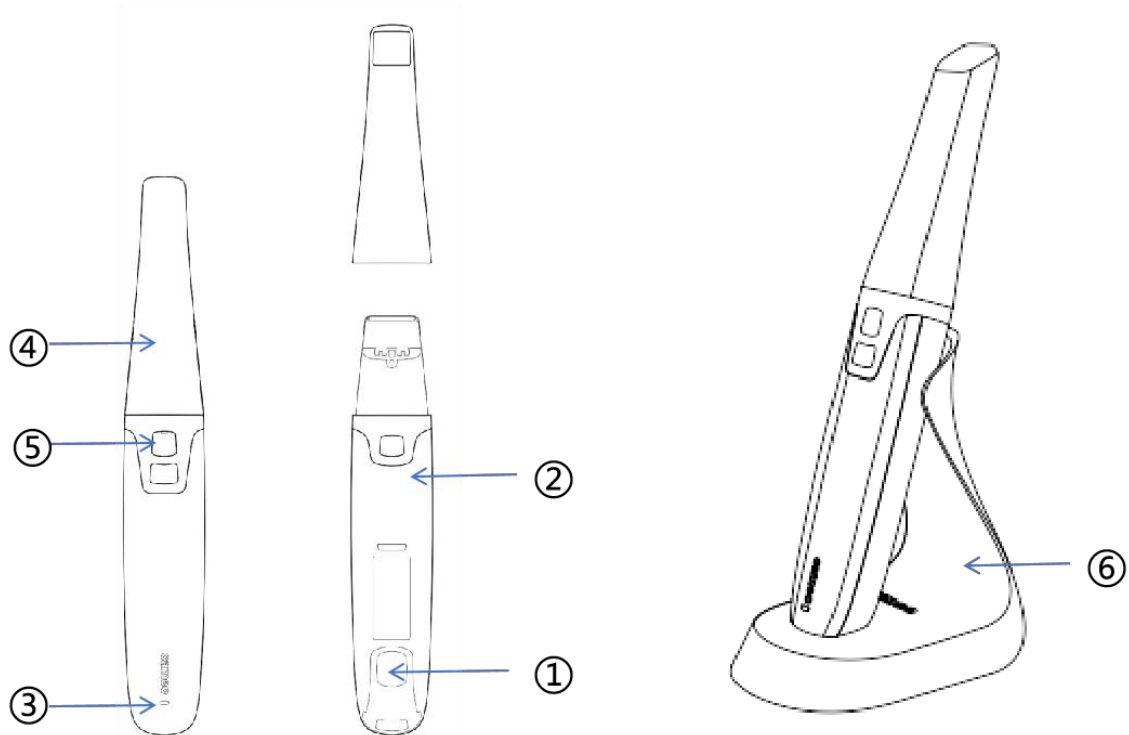

Step 1 : The lithium battery pack is connected to Handpiece Then the battery indicator starts to check the power level. The green light flashes for a few seconds and then turns into a steadily lit green light. Meanwhile, the status indicator's blue light flashes. Step 2: Plug in the scanner tip . When the connection to the wireless network scanner is successful, press the key 5. If the key indicator lights up green, it indicates that scanning can be carried out. Status of the key indicator: Normal connection/standby: The blue light stays on constantly; Scanning state: The green light stays on constantly; Abnormal connection: The blue light flashes; Sleep state : The cyan light stays on constantly. Step 3: After the scanning is completed, the scanner can be placed into the cradle, and the scanner will enter the standby mode and stop the scanning operation. Step 4: Disconnect the lithium battery pack from the scanner to cut off the power supply of the whole device.