

TPC 3D



Scanner Connecting of Infinity 5

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

Follow the steps below to complete the wired connection.

1. Remove the dusk cap



2. Install the scanner tip

Push the scanner tip firmly onto the scanner handpiece to ensure a secure connection.



3. Install the scanner charging base and place the scanner on the charging base.

Blue light flashing indicates that the scanner is waiting to connect.






4. Install the Wi-Fi adapter on to the computer and insert it in to the USB3.0 interface.




Attention: The Wi-Fi adapter shall support TP-Link protocols and the 5 GHz band, and hold an ANATEL certificate. Please download the TP-Link driver first when using it on a new computer. The address to download TP-Link driver: <https://oss.aiditecloud.com/upload/202404/file/1768366699188.zip>

5. Connect the wireless network interface controller in two ways:

□ Open the WiFi settings, to check if there are two wifi channels. If the WiFi settings do not display two WiFi channels, please install the driver and check whether the indicator light on the receiver is yellow.


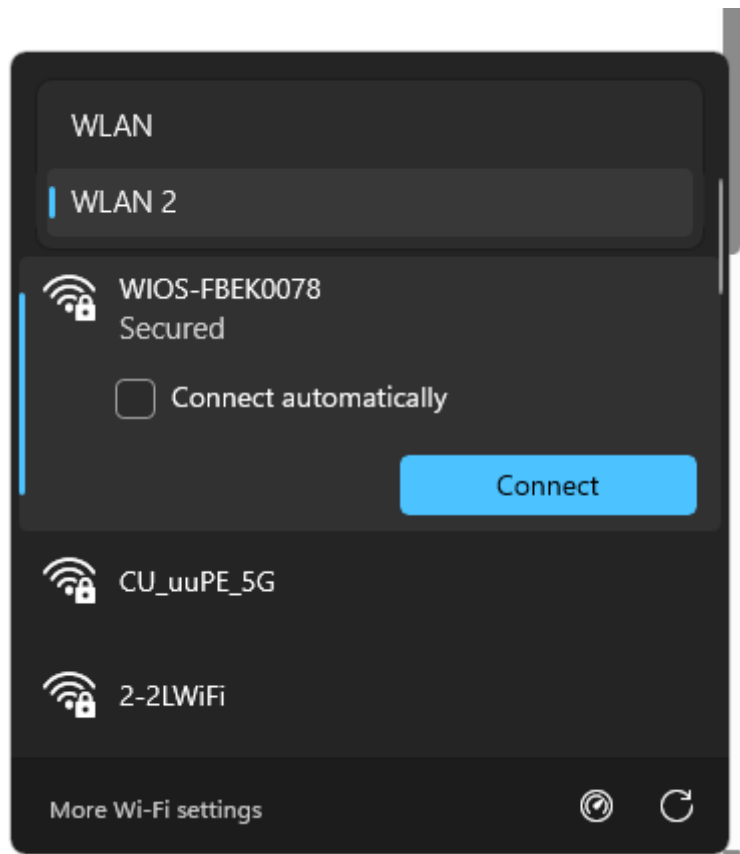


If no scanner signal beginning with "WIOS" can be found in the WiFi settings, you may need to reset the channel. • First, remove and reinsert the scanner battery. • Then press and hold both buttons on the front of the scanner. Wait until the indicator light on the scanner buttons changes from a blue breathing effect to a solid bright blue light . — this confirms that the channel reset is successful.

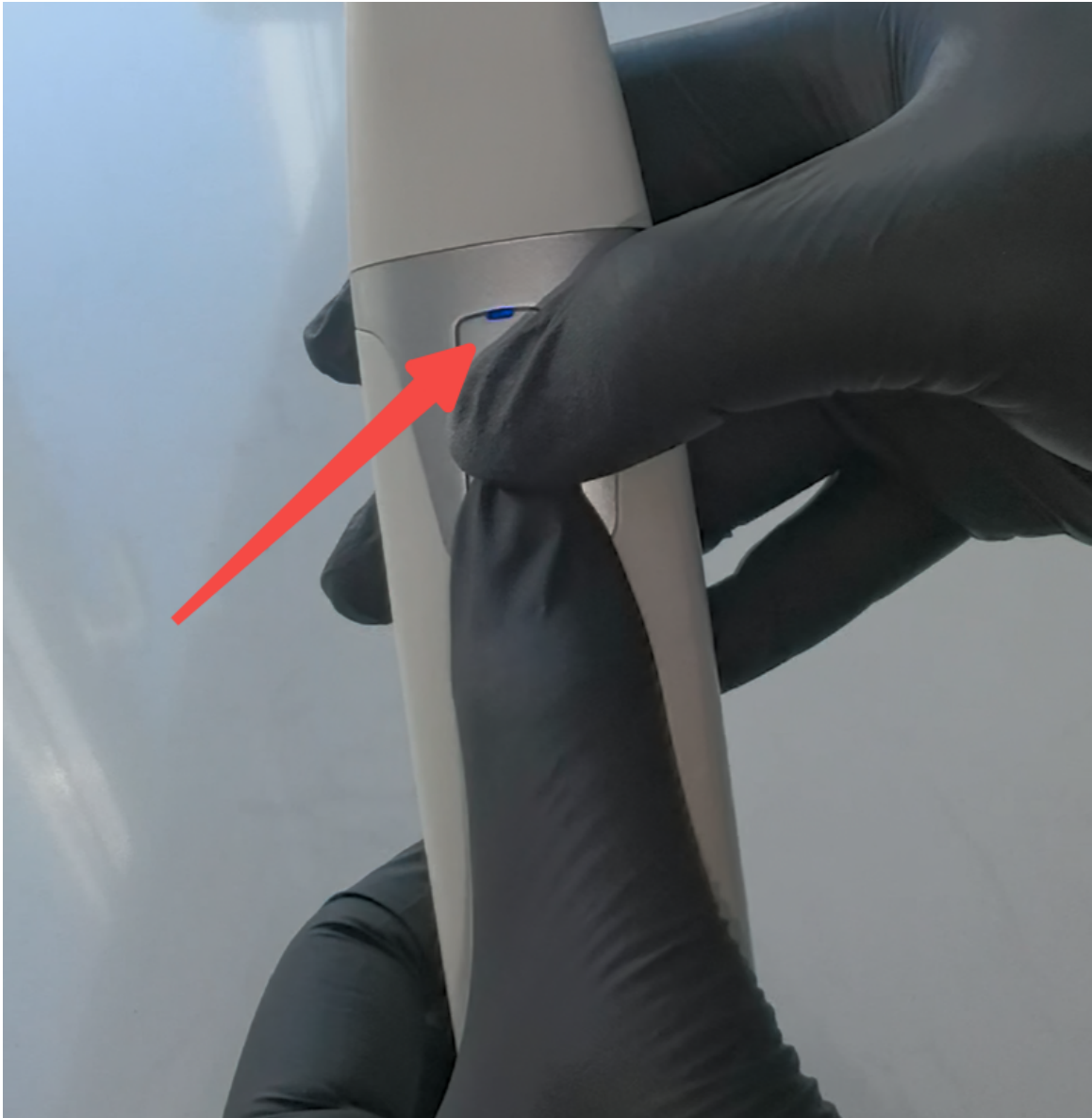



- Last, remove and reinsert the battery once more and place the scanner on the charging dock. Recheck the WiFi settings for the scanner signal. □ Open the software settings, click "Wireless device connection assistant", select the serial number of the corresponding device, and connect

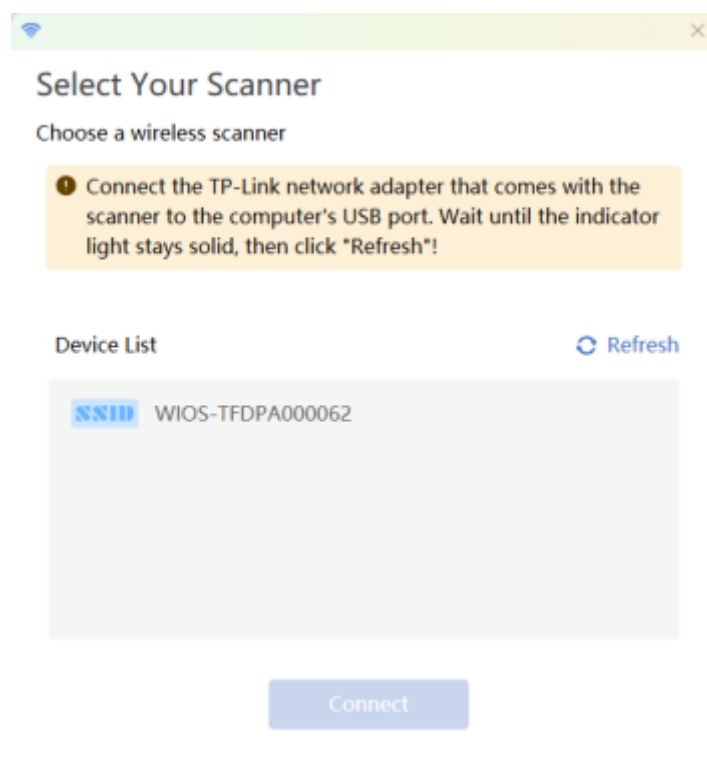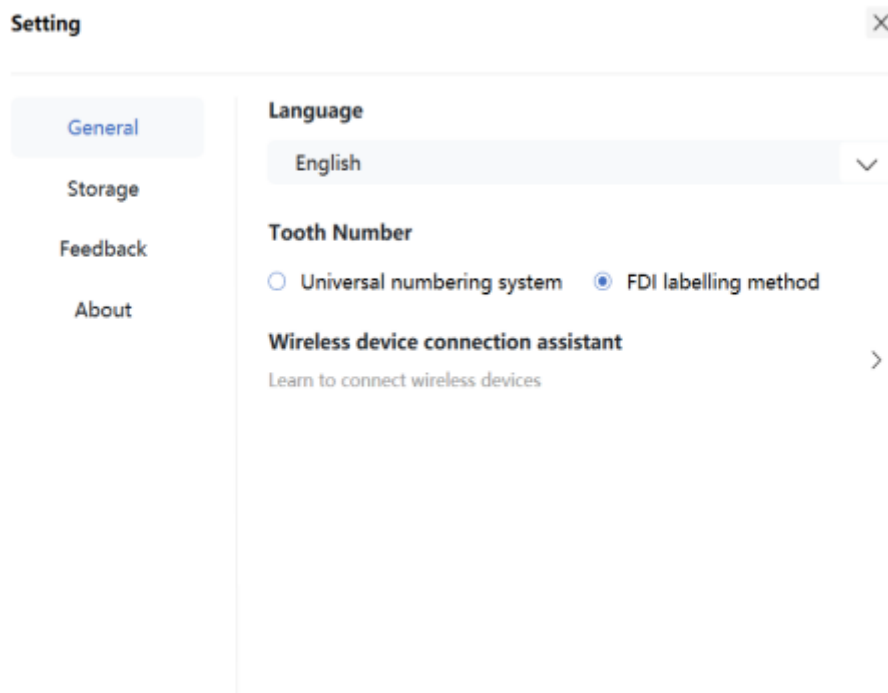

The scanner indicator is constantly on in blue, indicating a successful connection.