

TPC 3D



Add Patient and Case

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

2.1.1 Add patient

The screenshot shows the Aidite application interface. On the left, there is a sidebar with the 'Add Patient' button highlighted by a red box. A red line connects this button to the 'New Patient' form that is open in the center. The form contains the following fields:

- First Name ***: Text input with 'Test' entered.
- Last Name ***: Text input with 'Aidite' entered.
- Gender ***: Radio buttons for Male (selected), Female, and Others.
- Birthday ***: Date picker showing '3/15/2007'.
- Phone**: Text input.
- Email**: Text input.
- Buttons**: 'Cancel' and 'Confirm' buttons at the bottom.

The background interface includes a top navigation bar with 'All Cases', 'Scanned cases', and 'Empty Cases' tabs, and a '+ Add Cases' button on the right. A search bar and filter icons are on the left. At the bottom, there is a status bar with 'Total 0', '9/Page', and pagination controls.

2.1.2 Add case

- Select patient and add a new case.

The screenshot shows the Aidite application interface. On the left, there is a sidebar with the 'Add Patient' button. A list of patients is displayed below it:

- Aidite Test**: Age:36, 8/25/2026 13:50
- Test Aidite**: Age:19, 8/25/2026 13:47

In the center, the 'Select Patients' form is open. It contains a search bar and a 'Near Term Selection' section with two radio buttons: 'Aidite Test' (selected) and 'Test Aidite'. At the bottom are 'Cancel' and 'Confirm' buttons. A red line connects the 'Add Case' button (highlighted with a red box in the top right) to the 'Select Patients' form. The background interface is the same as in the previous screenshot, with the 'All Cases' tab selected.

- On the case creation screen, you can modify the following information.

Aidite

Test Aidite
Create Time 8/25/2026 13:52 **①**

Doctor * **②**
Test_doctor

Operator **③**
Test_operator

Due Date **④**
8/25/2026

Case Type

Note **⑤**
Please enter content:
Design & Processing Requirements:
Precautions:

⑥ (Tooth diagram)

⑦ (Tooth 15)

⑧ (Crown)
Pontic
Inlay
Veneer

⑨ (Tooth Shade)
VITA Classic Shade | A1

⑩ (Implant Type)
No Implant

⑪ (Materials)
Zirconia

⑫ (Preoperative Scan)
Yes **⑬**

Preview Export Send Exocad Scan

□ Doctor Information

Doctor *

Test_doctor

✓ Test_doctor

Doc1

Doc2

Doc3

Operator

Test_operator

✓ Test_operator

Operator1

Operator2

Operator3

☐Operatot Information ☐Due Date

Due Date

8/25/2026

◀ August 2026 ▶

S	M	T	W	T	F	S
26	27	28	29	30	31	1
2	3	4	5	6	7	8
9	10	11	12	13	14	15
16	17	18	19	20	21	22
23	24	25	26	27	28	29
30	31	1	2	3	4	5

Implant

Health Report

Case Type

Restoration

Orthodontics

Removables

Implant


Health Report

☐ Case Type ☐ Add Note ☐ Adult or Child






□Tooth Number You may select the tooth number corresponding to the working tooth. Furthermore, you may toggle between different tooth□numbering systems within the Settings menu.

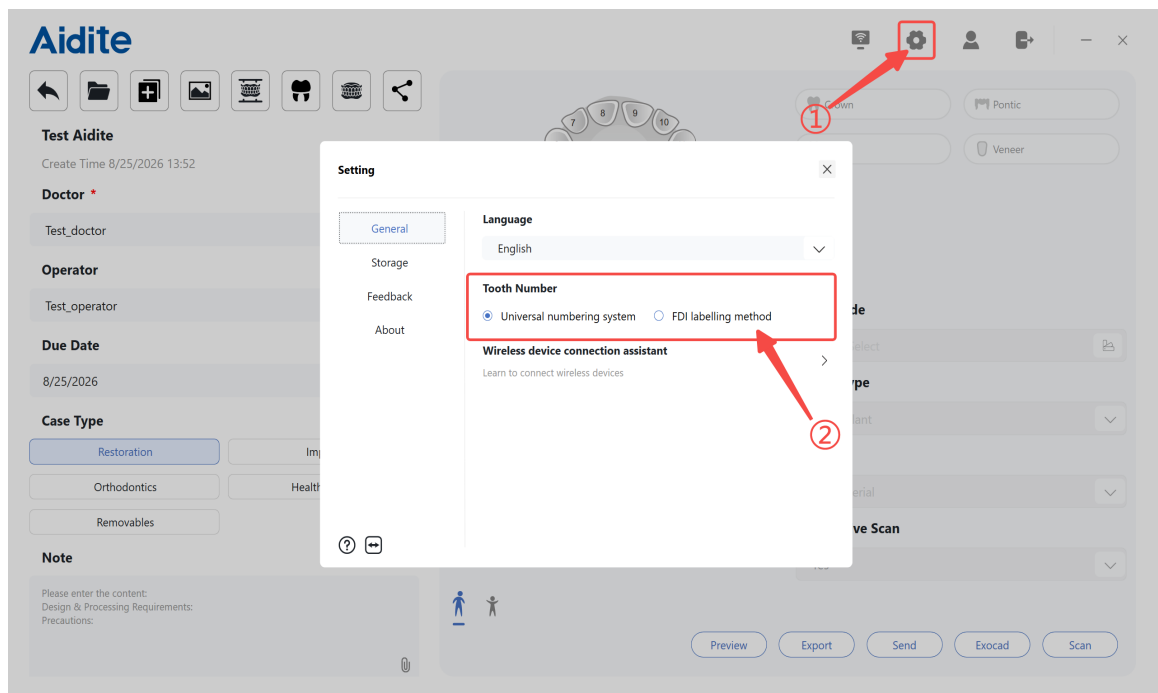

□Restoration/Implant/Orthodontics/Removalbe Type

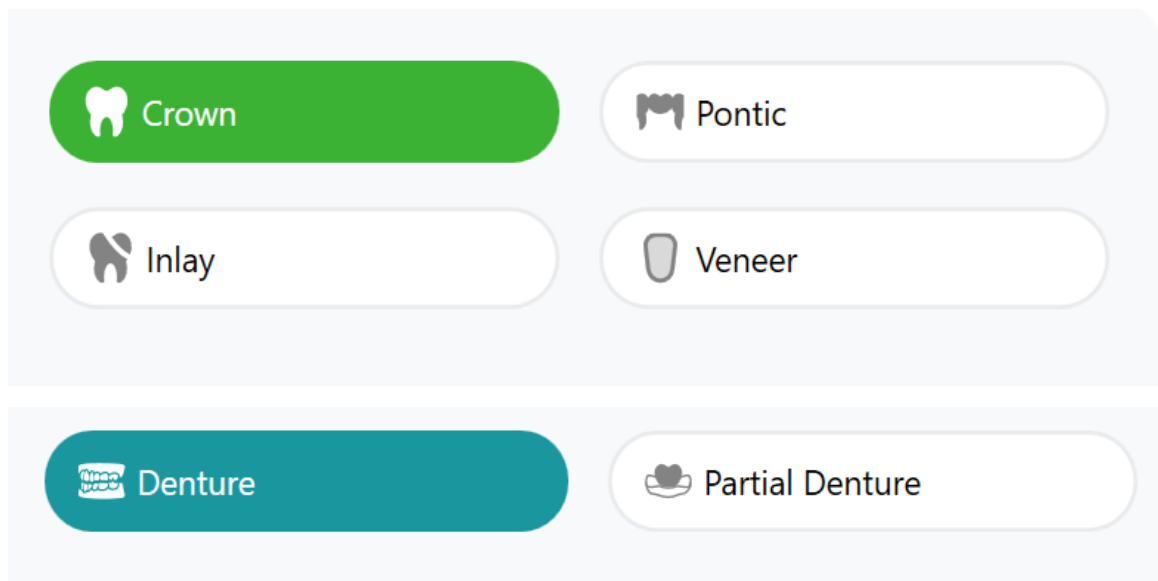

Implant Type

No Implant



✓ No Implant

Custom Abutment

Custom Abutment(manual)

Screw Retained

Screw Retained(manual)

□ Implant Type □ Materials

✓ Zirconia

GlassCeramic

Lithium Silicate

WAX

Titanium

Feldspar

NP Metal

PBC

Zirconia



Preoperative Scan

No



✓ No

Yes

□ Prepreative Scan or Not After export, you will get a PDF file of the case information. Example:

Case information report

Aidite

Clinic Name: aiditetest_lab2

Return Phone:

Return Address:

Doctor's Name: Test_doctor

Operator's Name: Test_operator

Due Date: 2026-08-25

Patient Name: Test Aidite

Birthday: 2007-03-15

Visit Date: 20260825154449

Scanning
process:

Restoration

Tooth Number: 15


Repair Type: Crown

Materials: Zirconia

Implant: No Implant

Tooth Shade: VITA Classic Shade | A1

Note:

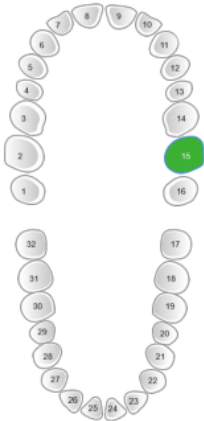

1