


TPC 3D



Introduction to the Buttons after Process

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

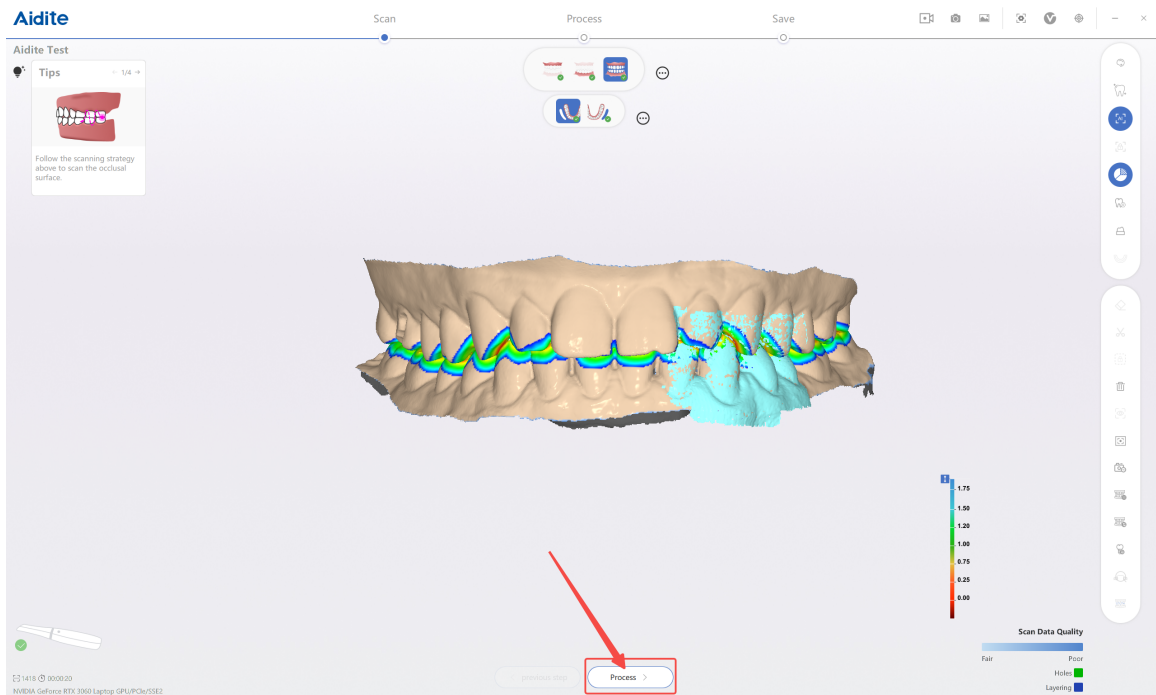

After clicking the process button, the toolbar will change.

Icon Name Function

True Color Turn on to display real colors of the scanned model.



Vivid Mode Turn this option on to display more vivid scanned images.



Smooth Enable this button to smooth the mesh.



Hole filling When enabled, the acquisition software will automatically fill holes in the scan. Eraser Move the eraser to the area you want to delete, and click to



delete the area covered by the green circle. Cut Move the cursor and press the left mouse button to draw a closed



curve. Release the mouse and click "OK" to delete the selected area. Reset View Restore to default size and position.



View HD Click this button to view HD image.



Image Undercut Click the button to inspect undercuts.



Check Occlusion Click this button to check the occlusion




Distance

True color scanning


Turn on to display real colors of the scanned model.




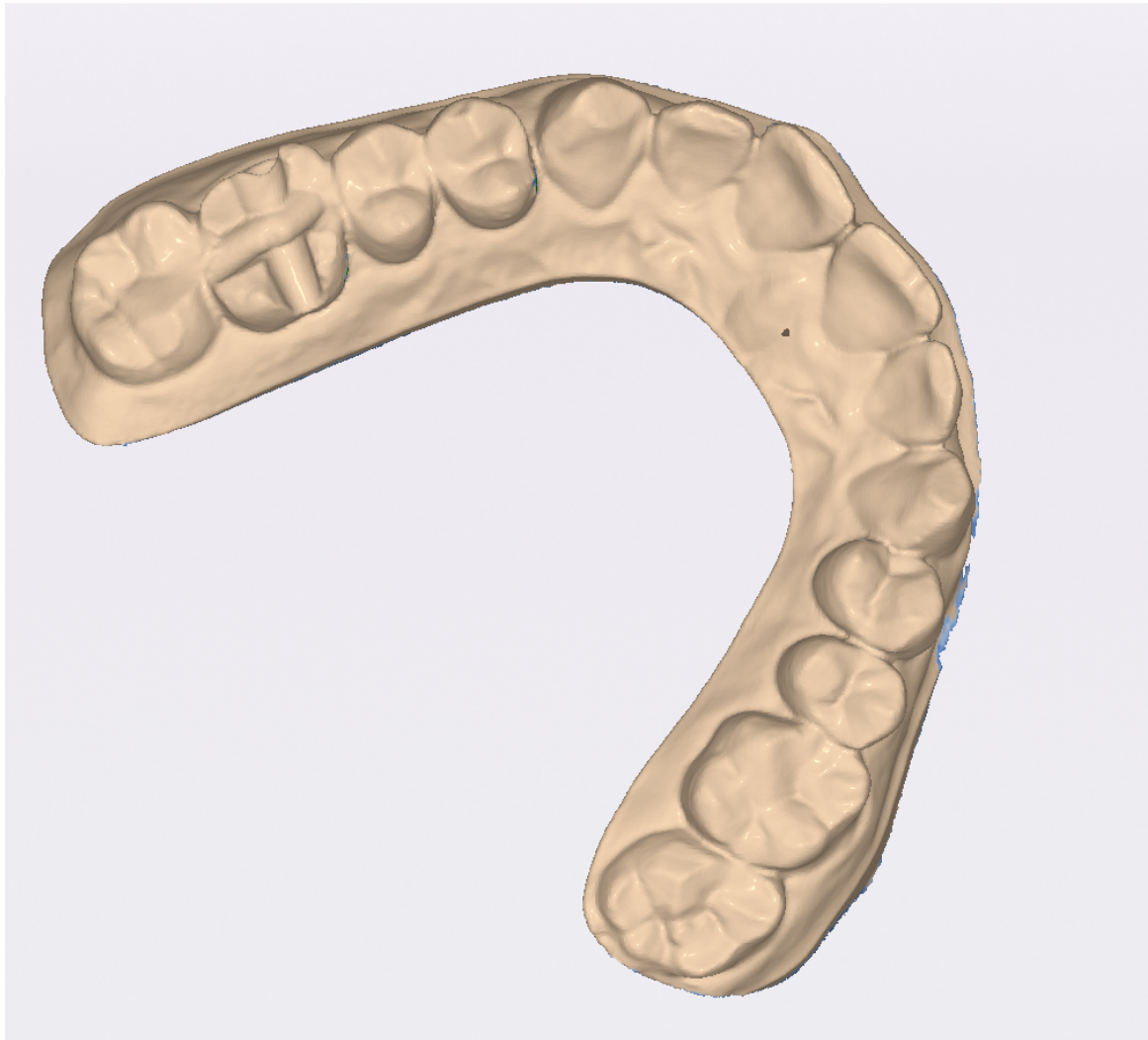

Vivid mode

Turn this option on to display more vivid scanned images.

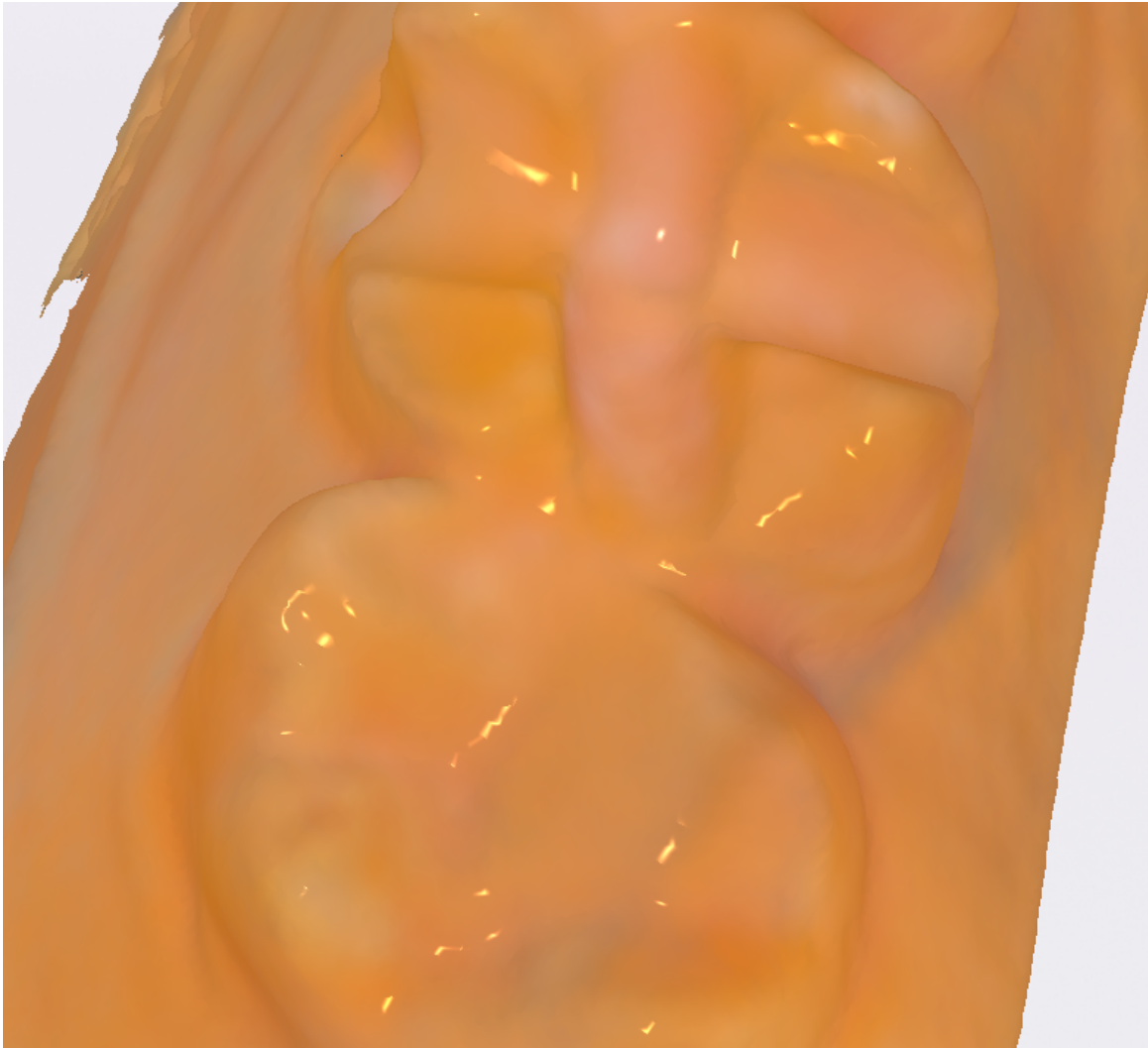

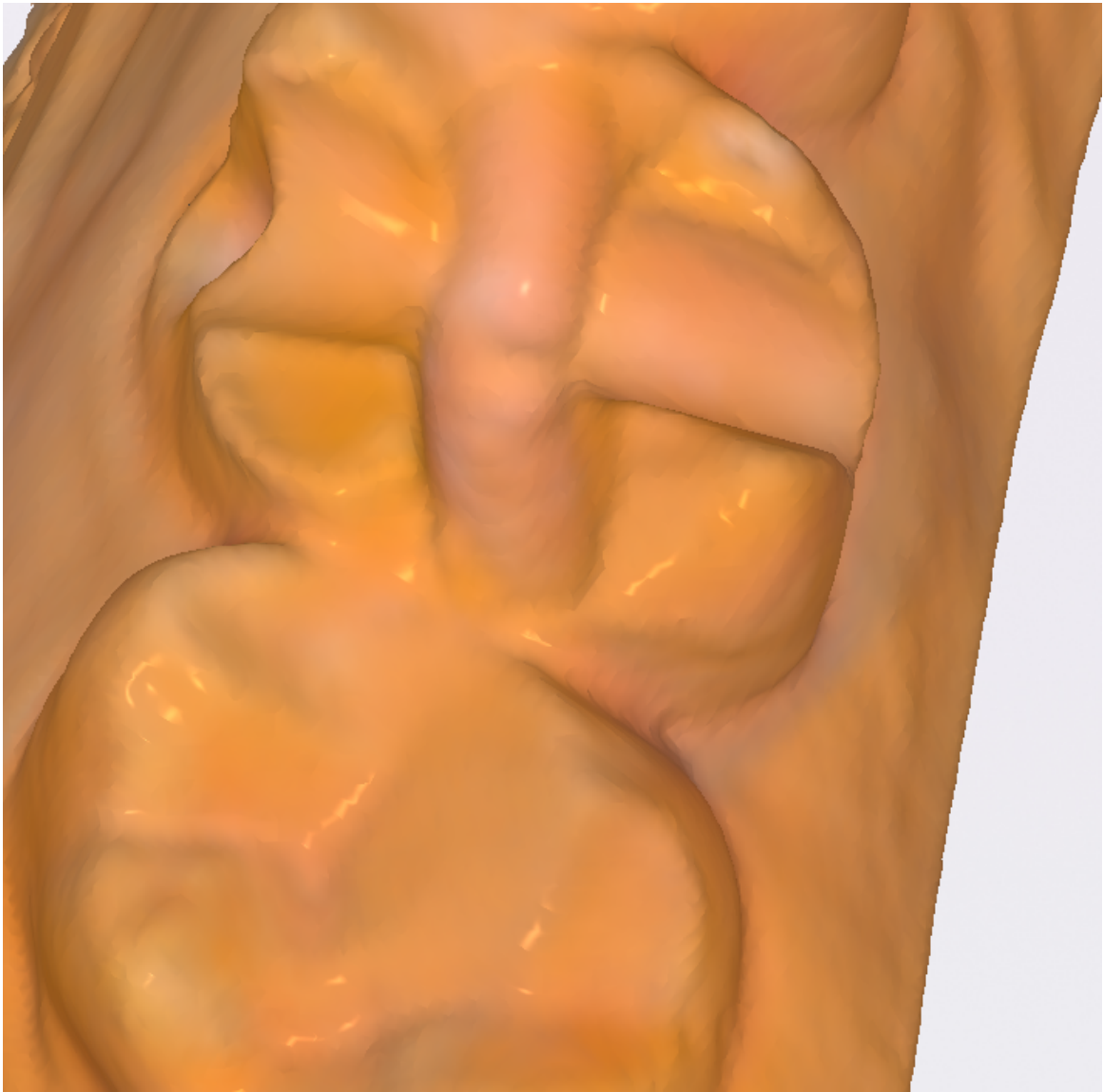



Smooth

Hole Filling

Eraser


Move the eraser to the area you want to delete, and click to delete the area covered by the green circle.

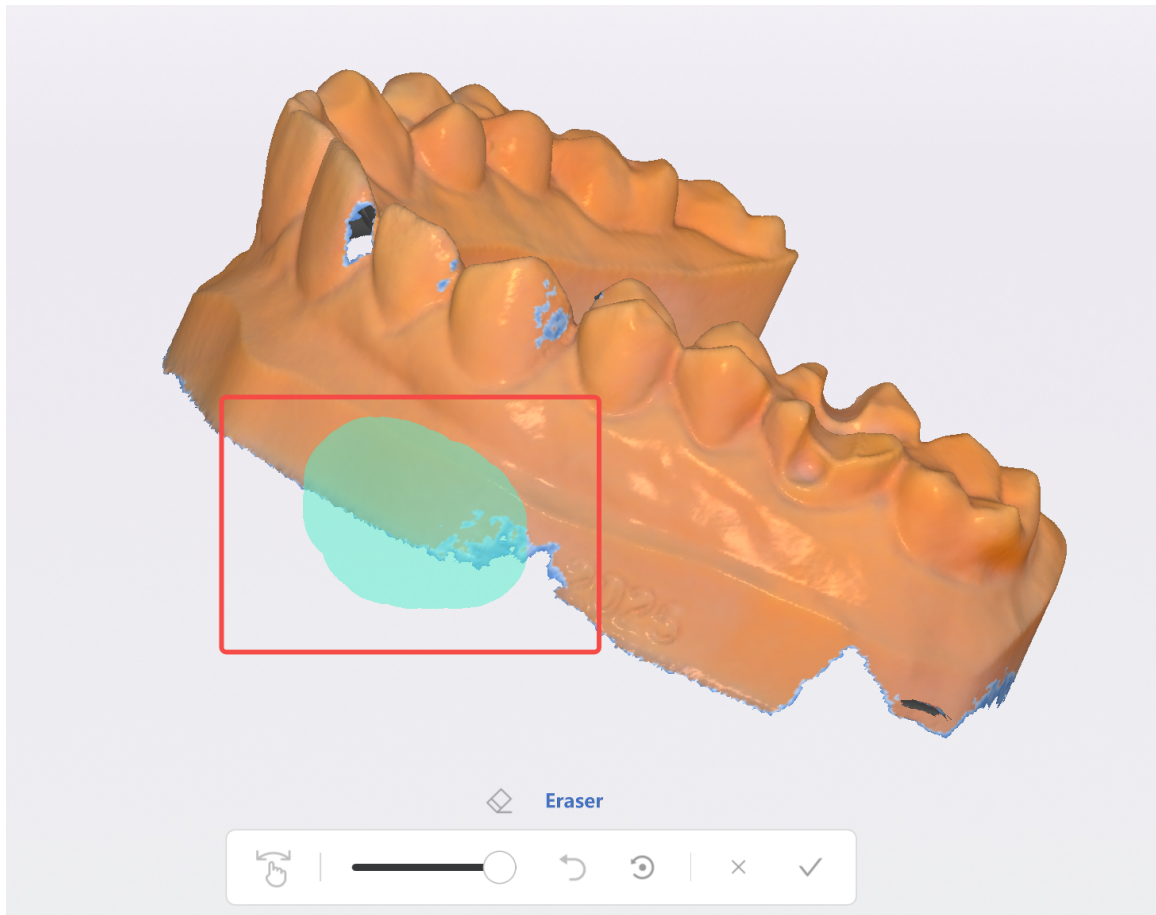


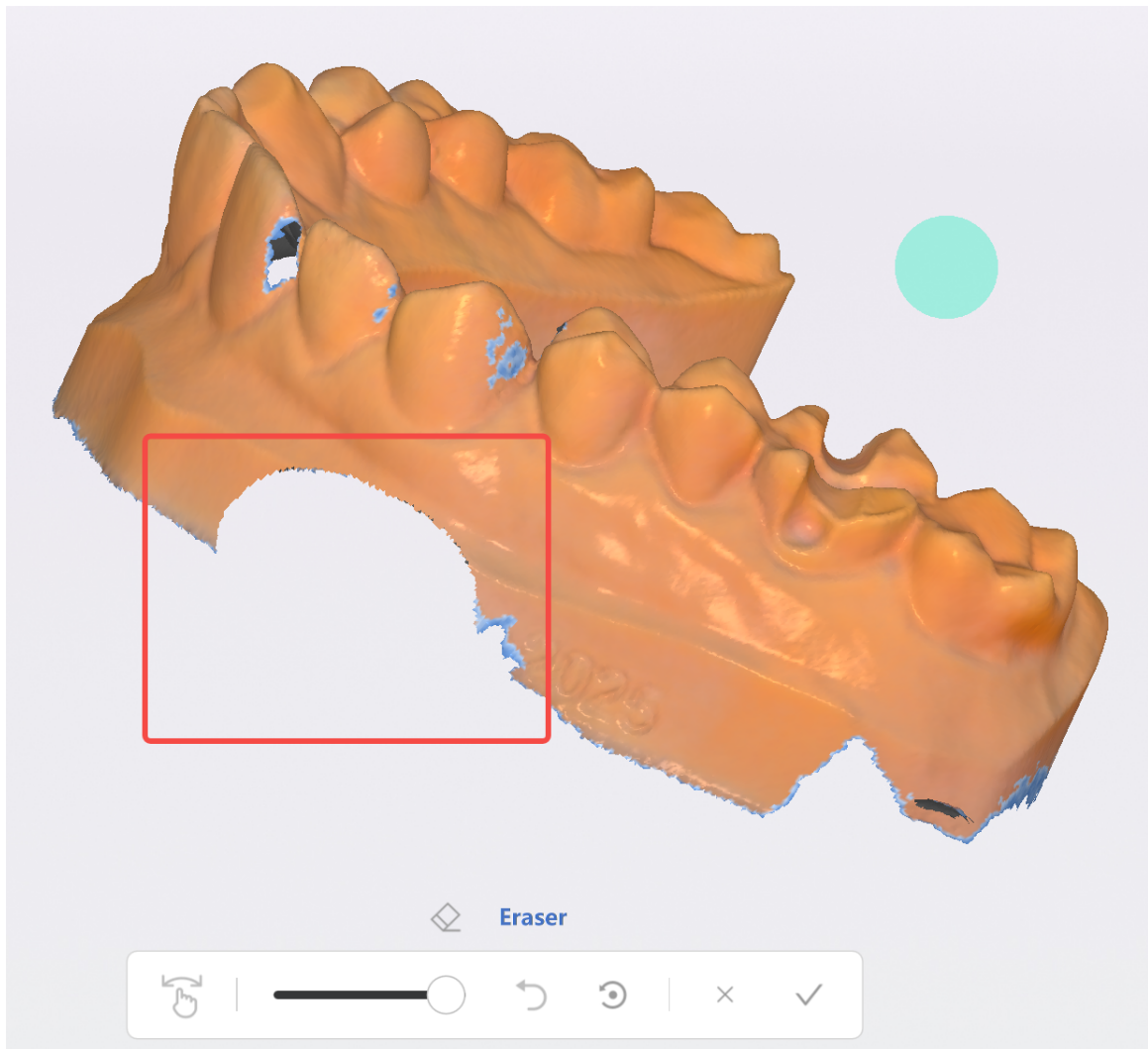

Cut

Move the cursor and press the left mouse button to draw a closed curve. Release the mouse and click "OK" to delete the selected area.

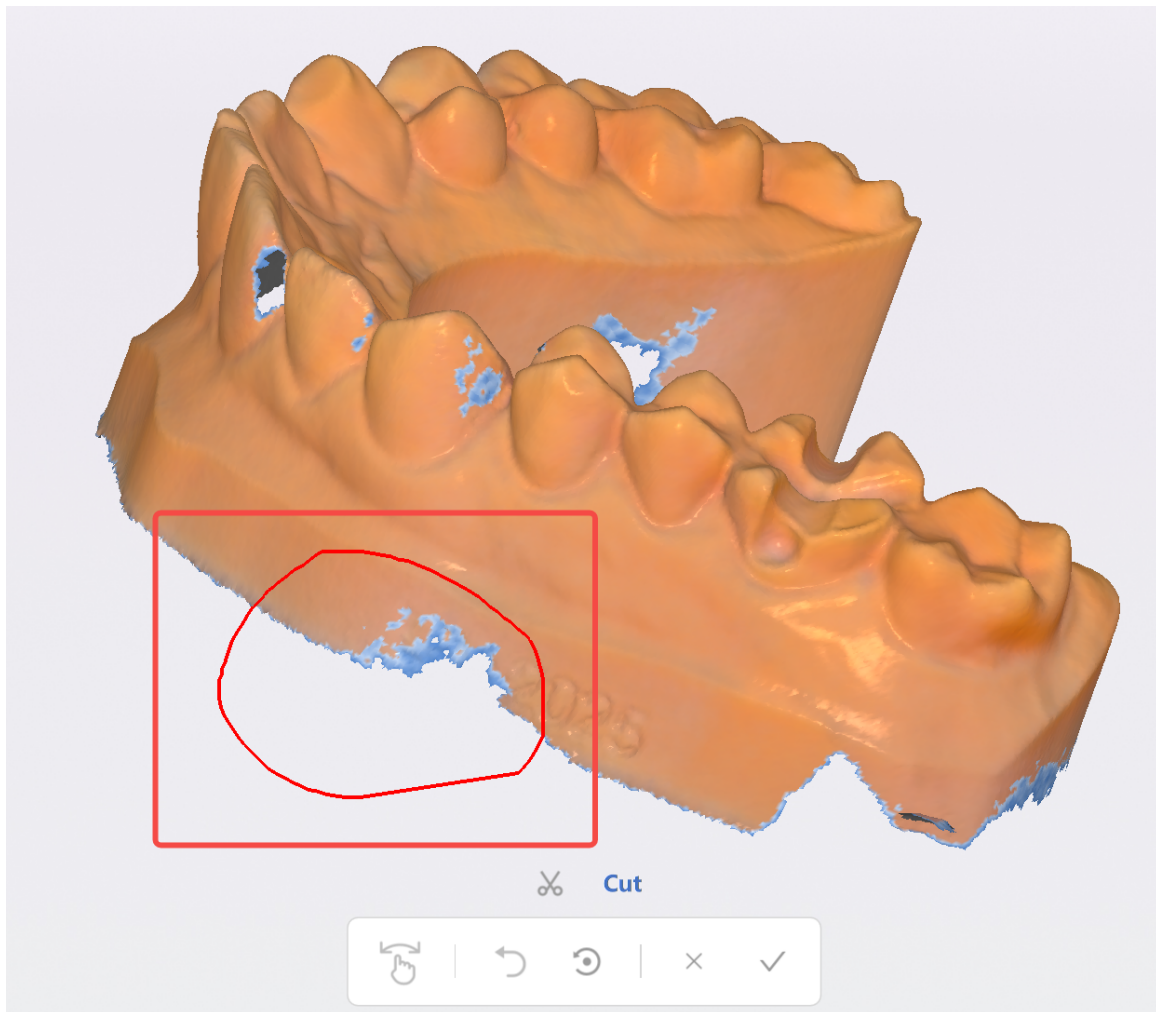


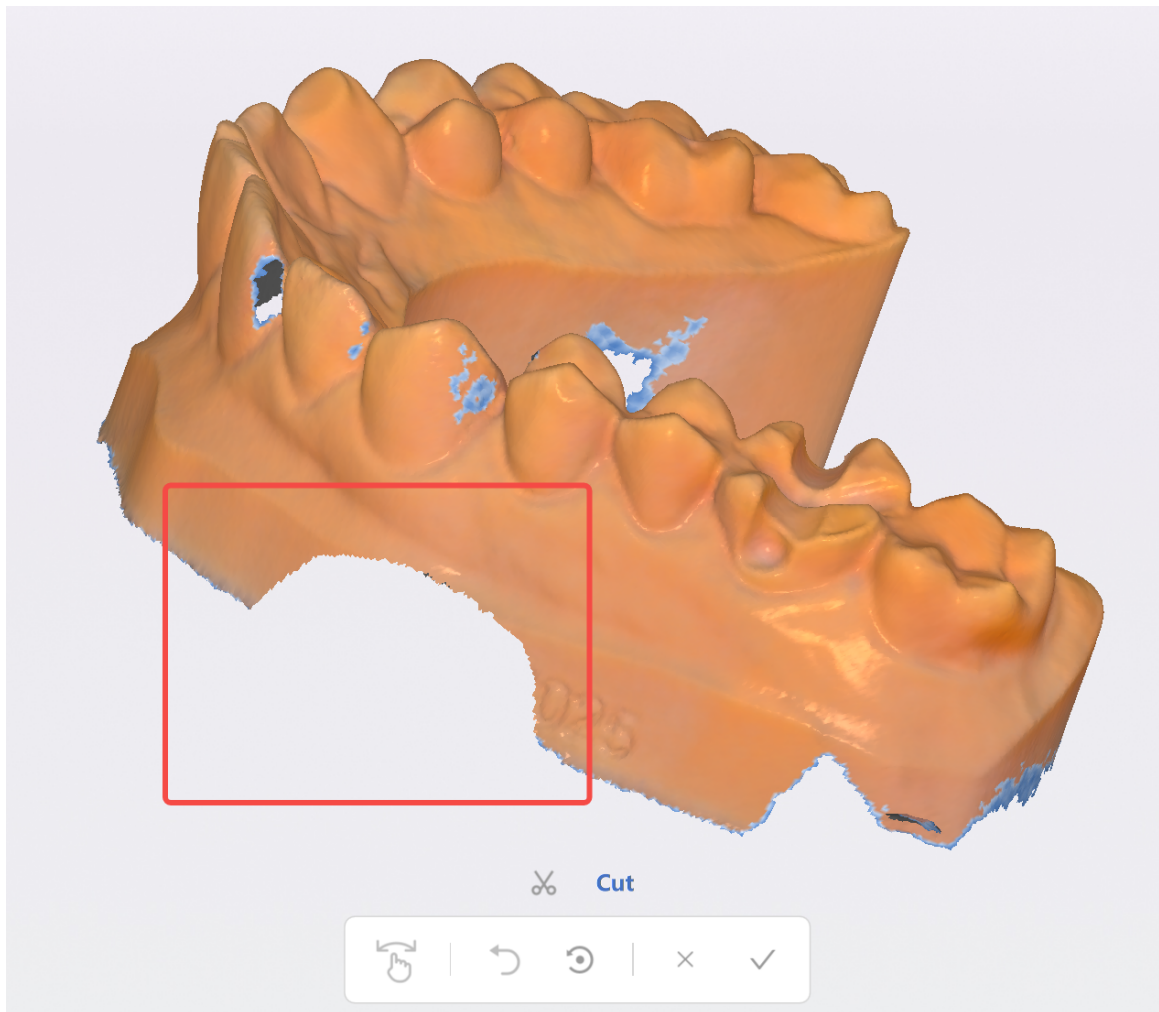


Reset View

Restore to default size and position.

Undercut Check

Click the button to inspect undercuts.

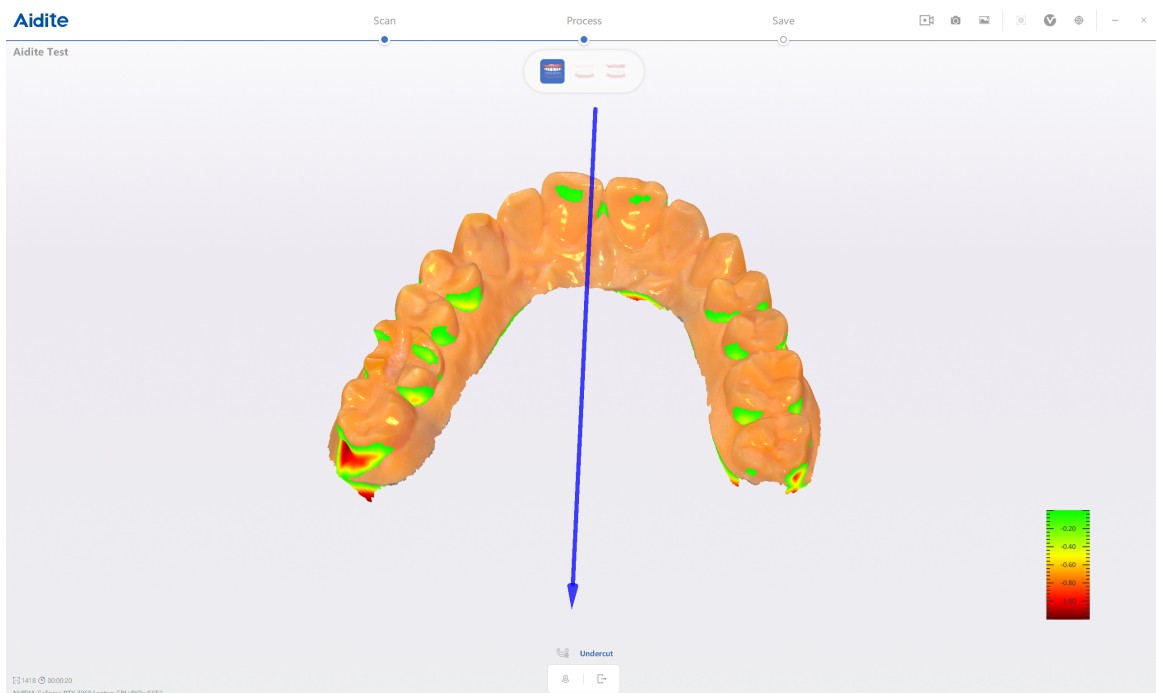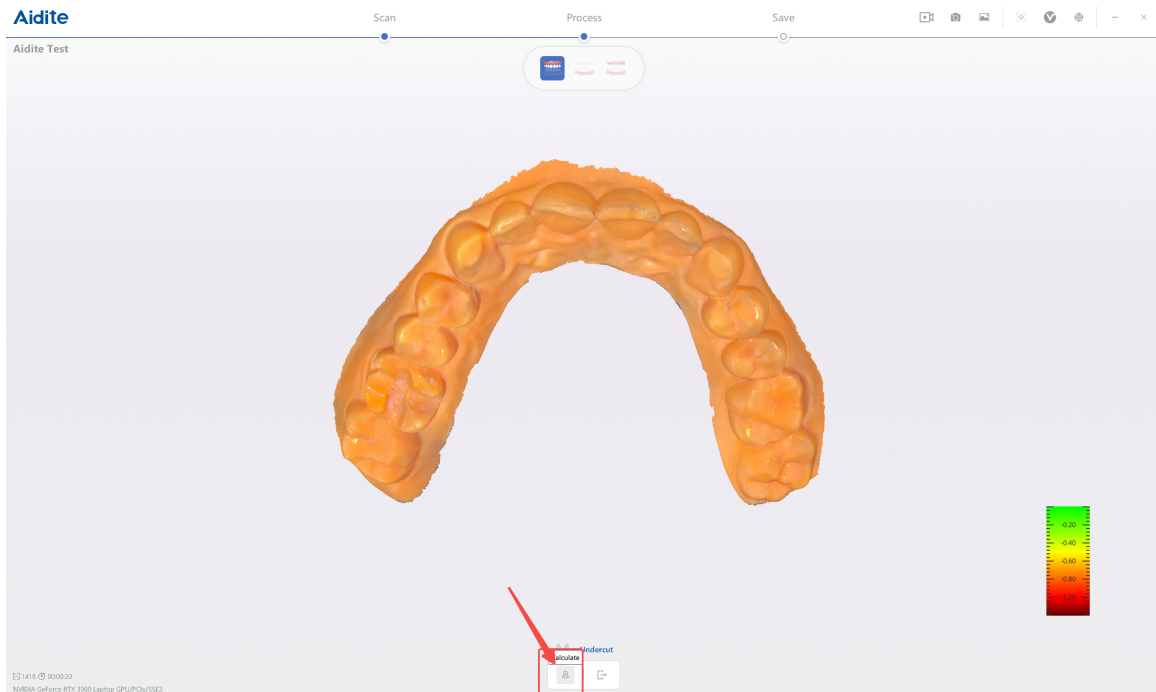

Occlusion Distance

Click this button to check the occlusion You can click split to split the upper jaw and lower jaw and drag the slider to check the occlusion

