

TPC 3D




Scanning Strategy

Intraoral Scanner User Guide


Everyone with a healthy and beautiful smile

This feature is designed to simplify the scanning process and help users capture more accurate scan data. The feature identifies scanning issues and provides appropriate guidance on how to revise them.


1. Follow the scanning strategy below to scan the molar area.

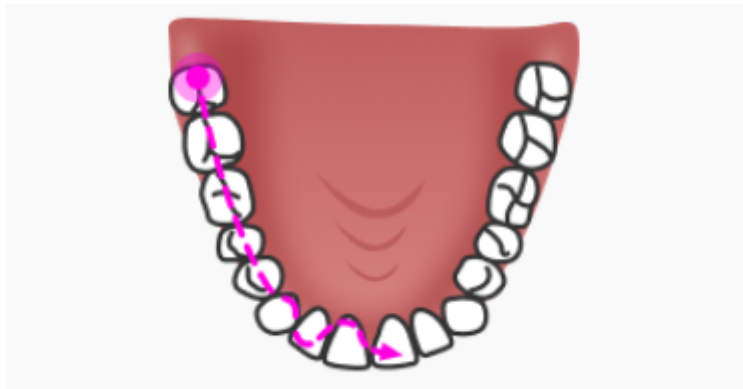


2. Maintain an appropriate scanning distance; do not stay too close to or too far from the teeth.

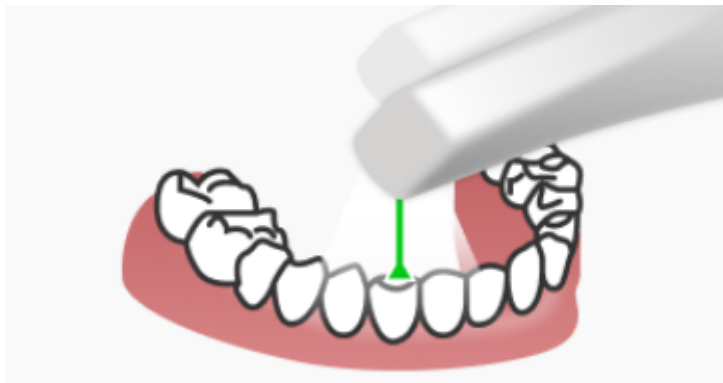


3. Ensure the scanning area is contiguous. If the scan is interrupted, return to the previously scanned area to resume.

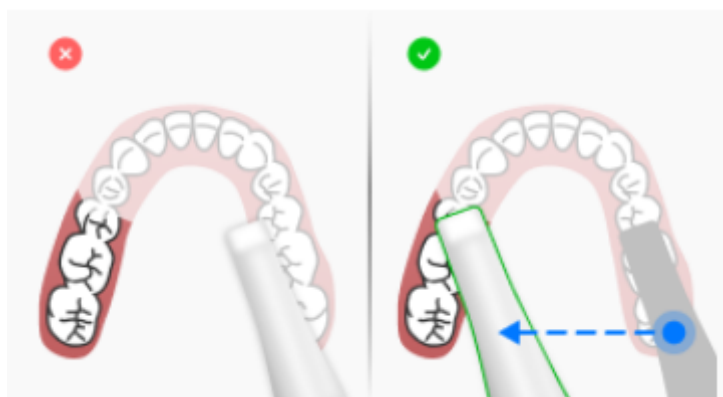

4. You can capture additional scan data, such as bridges, wax-ups, etc.

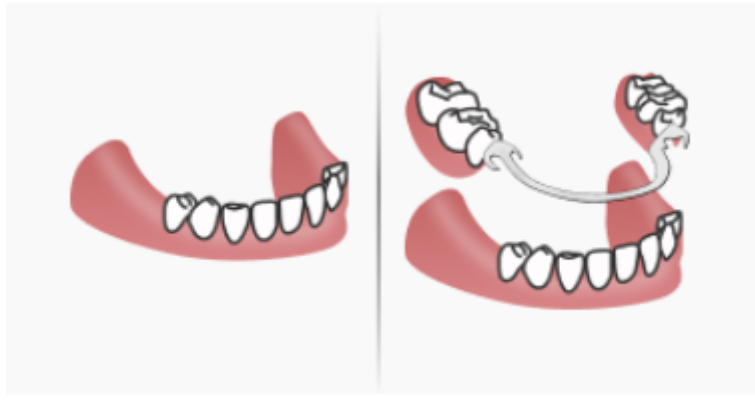

5. Please capture more data to fill the mesh gap. To ensure scan quality, you could lock the scanning angle.

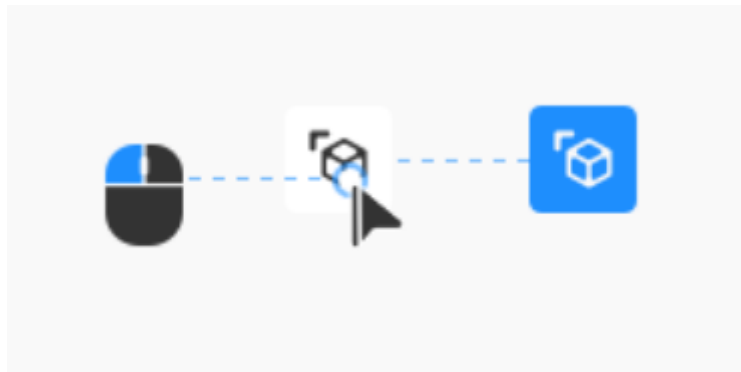

6. Please select a target area to mark. More refined optimization of the marked area will be made when clicking Next.

