

TPC 3D





Removable Cases Scanning Procedure

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

1. Add a patient and case

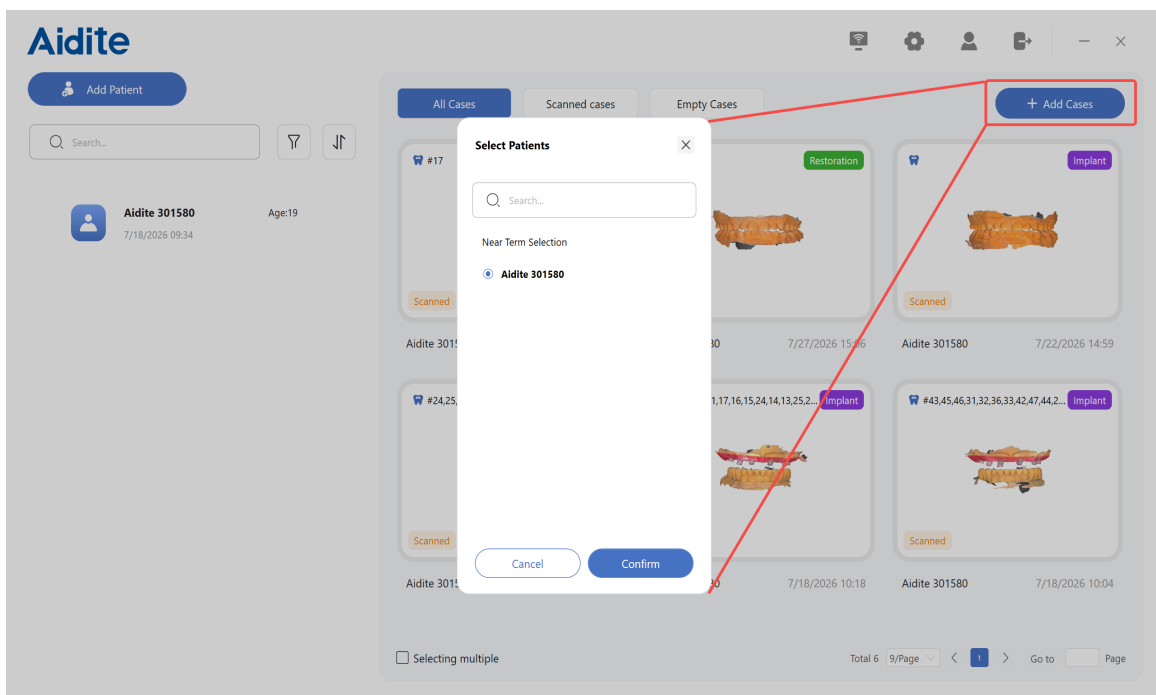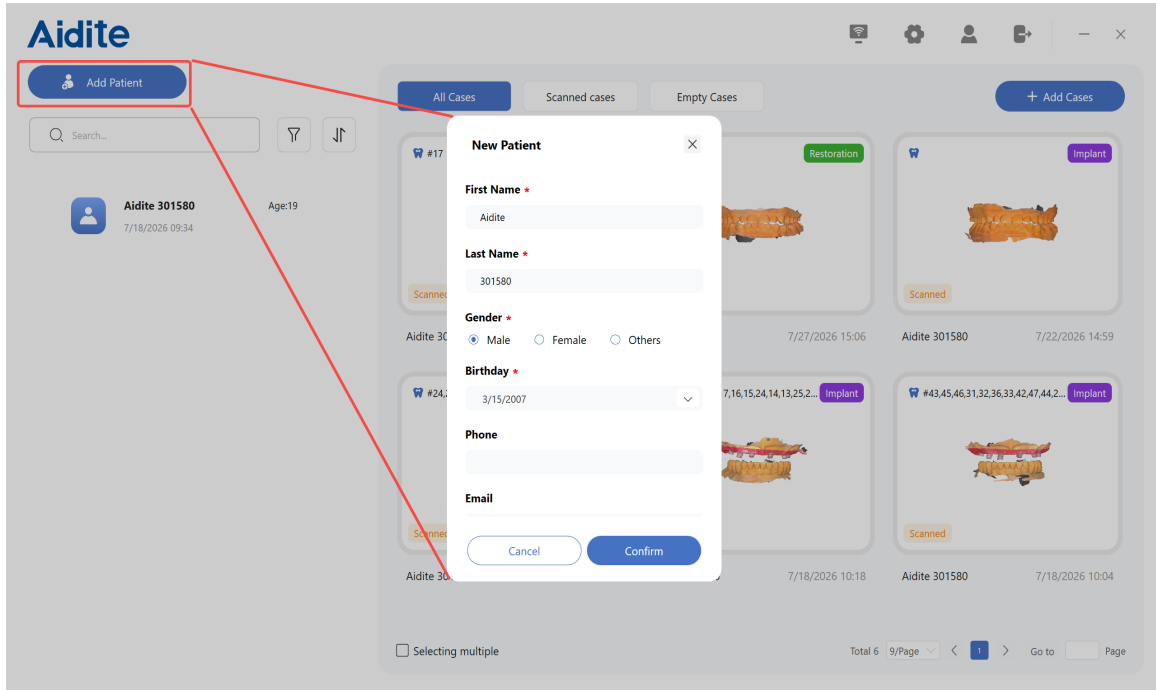

2. Set up case information.

Aidite

云平台 测试
Create Time 8/8/2026 13:27

Doctor *
111

Operator
Please Select

Due Date
Please Select

Case Type

Restoration Implant
Orthodontics Health Report
Removables

Note
Have a good day!

Denture Partial Denture

Upper
Lower


Method
Full Denture

Materials
No Material


Tooth Shade
Please Select

Preview Export Send Exocad Scan


3. Scanning:

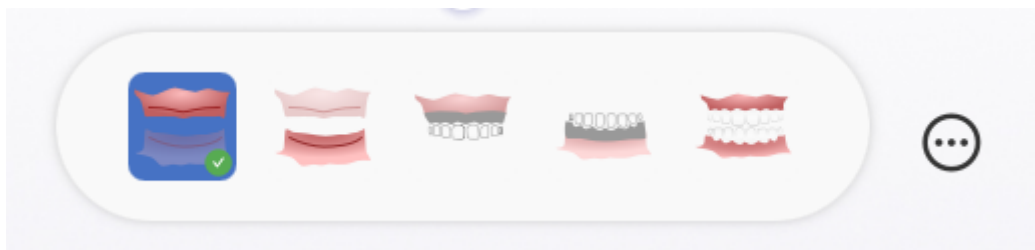


First, scan the edentulous jaw.

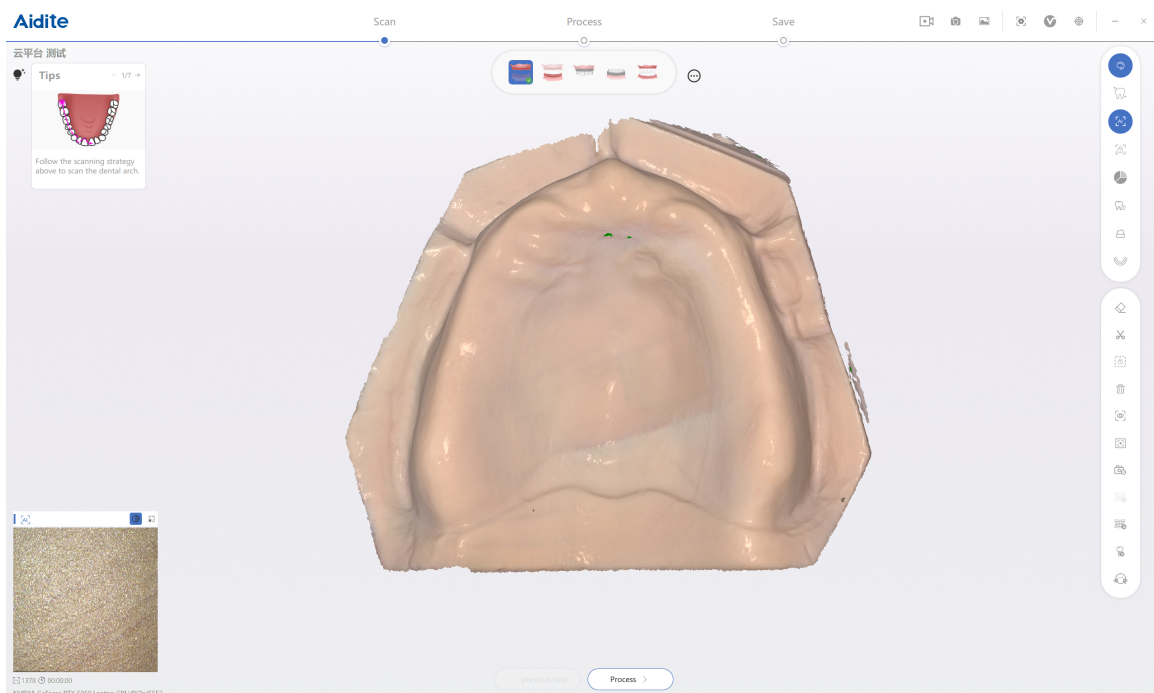


Second, scan the old denture.

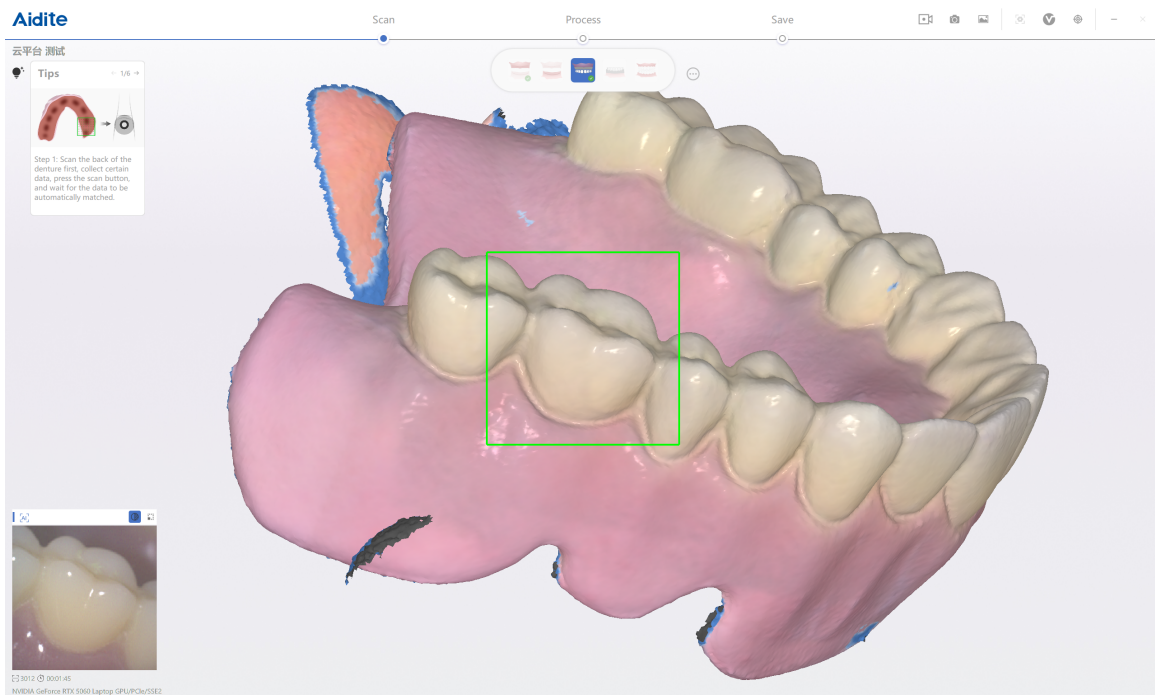

Then, the scanning software will automatically align the existing denture with the soft tissue of the edentulous arch.

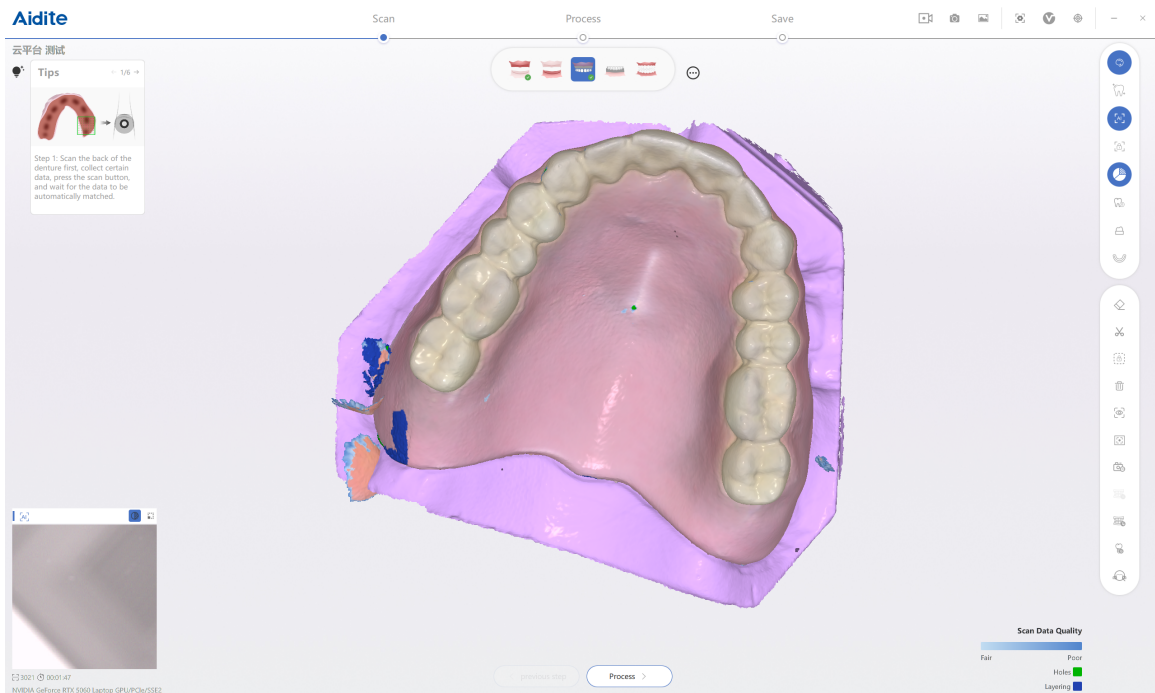

Third, scan the opposing arch.

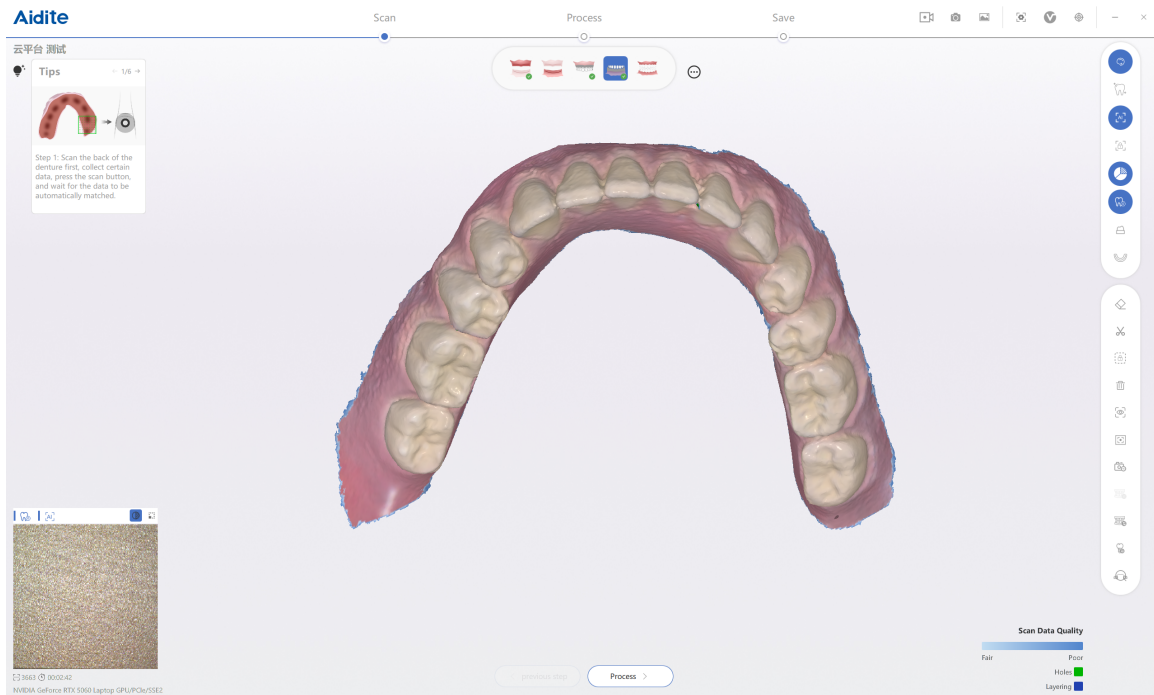

Last, scan the occlusion.

