

TPC 3D



Orthodontics Cases Scanning Procedure

Intraoral Scanner User Guide

Everyone with a healthy and beautiful smile

1. Add a patient and case

The screenshot shows the Aidite software interface. On the left, there is a sidebar with the 'Add Patient' button highlighted by a red box. A red line connects this button to the 'New Patient' form that is open in the center. The form contains the following fields:

- First Name ***: A text input field with the value 'Aidite'.
- Last Name ***: A text input field with the value '301580'.
- Gender ***: Radio buttons for 'Male' (selected), 'Female', and 'Others'.
- Birthday ***: A date picker showing '3/15/2007'.
- Phone**: A text input field.
- Email**: A text input field.

At the bottom of the form are 'Cancel' and 'Confirm' buttons. The background shows a grid of dental cases with labels like 'Restoration', 'Implant', and 'Scanned'.

The screenshot shows the Aidite software interface. On the right, there is a '+ Add Cases' button highlighted by a red box. A red line connects this button to the 'Select Patients' form that is open in the center. The form contains the following elements:

- Search...**: A search input field.
- Near Term Selection**: A list of patients, with 'Aidite 301580' selected and indicated by a blue dot.

At the bottom of the form are 'Cancel' and 'Confirm' buttons. The background shows the same grid of dental cases as the first screenshot.

2. Set up case information.

Aidite

云平台 测试

Create Time 8/13/2026 07:50

Doctor *

111

▼

✎

Operator

Please Select

▼

✎

Due Date

Please Select

▼

Case Type

Restoration

Implant

Orthodontics

Health Report

Removables

Note

Please enter the content:

Design & Processing Requirements:

Precautions:

Fixed Orthodontics

Invisible Orthodontics

Removable Orthodontics

Orthodontics Stages

Startup

▼

Preview

Export

Send


Scan

3. Scanning:


Orthodontics Cases Scanning Procedure

2

First scan the upper jaw

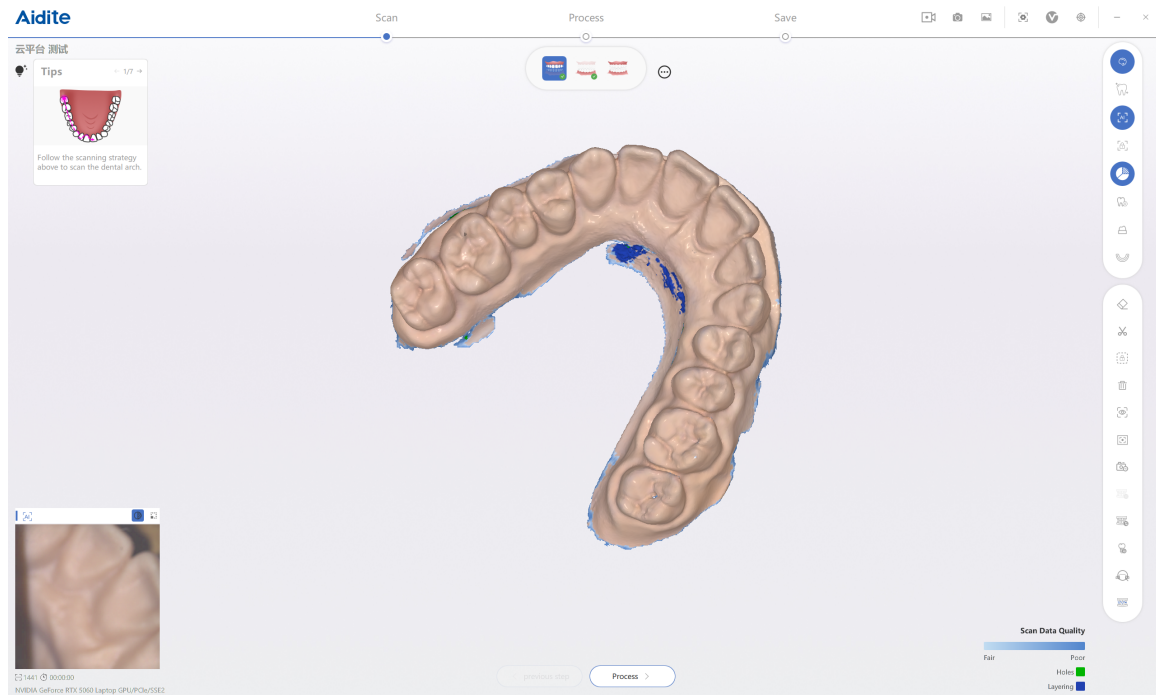

Second scan the lower jaw

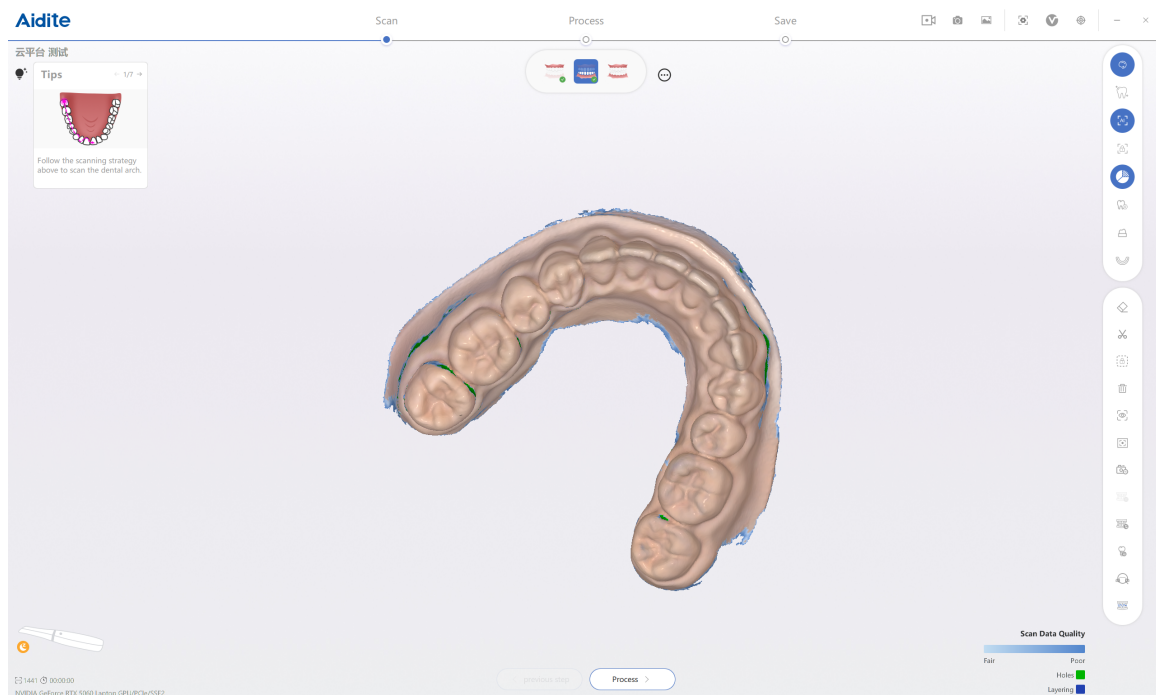

Finally scan the occlusion

