

# TPC

# Mirage

Chair Mounted Delivery System



Model #2000




### Packing Dimensions

Unit:	51" x 21" x 17"	60 lbs.
PMU:	23" x 19" x 21"	36 lbs.
Cuspidor:	18" x 14" x 14"	20 lbs.



### Cuspidor Bowl:

The deep contoured bowl is made of fine porcelain for easy cleaning and infection control. The cuspidor operates smoothly with the timed bowl rinse, cup filler and swing-away arm.

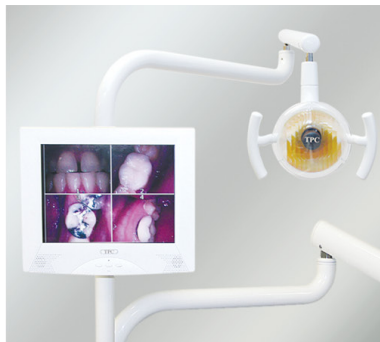

### Optional:

Multi media flat screen monitor for chair side intraoral imaging, patient education and DVD viewing.



### Control Head:

Easy to read pressure gauge. Easy access controls for simple adjustments. Ultra soft handpiece tubing minimizing hand and wrist stress. 3 handpiece automatic air diaphragm block

## Standard Features

- Asepsis automatic control for three handpieces
- Individual handpiece pressure adjustments
- Non-retracting water coolant with on/off
- Individual water coolant flow controls for all handpieces
- Air coolant flow control
- Handpiece tubing flush system
- Ultra soft handpiece tubing
- 4-hole straight gray asepsis tubing with Midwest connector
- Stainless steel tray with non-skid pad
- Disc foot control wet/dry toggle with 7' of smooth tubing
- Equipment utility center with asepsis housing
- Vacuum canister with solids collector screen
- Water outlet quick disconnect with volume control
- Air outlet quick disconnect for prophy jet
- Self contained water system with city/bottle selector switch
- Straight light post for unit mount dental light
- Pre-wiring for fiber optic light source
- Air-actuated counterbalanced arm lock
- 3 Way Syringe on unit head

## Gravity Drain Cuspidor

- Contoured bowl is made of fine porcelain for easy cleaning and infection control. The cuspidor operates with adjustable timed bowl rinse and cup filler
- Swing away cuspidor bowl with safety breakaway feature

## Assistant's Instrumentation on Breakaway Telescoping Arm

- 1 autoclavable HVE with asepsis tubing
- 1 autoclavable saliva ejector tubing
- Breakaway feature prevents damage if the telescoping arm contacts the assistant's stool as chair is being lowered

## Options

- Mounting bracket to mount flat screen monitor on light post
- Super thin X-ray viewer to mount on unit control head
- Built-in Led Curing Light
- Built-in Intraoral Camera
- Built-in Scaler
- Illumination Light Source
- Prophy Jet
- Post-mounted plastic utility tray
- Unit head mounted touch pad control

# Mirage

## Chair Mounted Delivery System

**TPC**  
**5 Year**  
**Limited**  
**Warranty**



Model #2015

## Packing Dimensions

Unit:	51" x 21" x 17"	60 lbs.
PMU:	23" x 19" x 21"	36 lbs.



Model #2000



### Standard:

Built in self contained water bottle system including city/bottle selector switch.



### Standard:

Built in air and water quick disconnect for ultrasonic scaling and prophy jet.