

It's All About Comfort.

- Thin backrest
- Double articulating headrest
- Easy to clean ergonomic design
- Easily programmable foot control with 3-position memory
- Cast aluminum construction for strength and durability
- Lifting capacity: 650 lbs.



- Chair Rotates 60°
(30° each side of center)



- 14" low base position enhances operator comfort for semi-supine procedures



- Functional armrests rotate to permit easy entry/exit for wheelchair patients

Standard Features:

- Thin back with lumbar support
- Anatomically correct back pivot point
- Adjustable double articulating headrest with contoured cushion
- Four manual chair movements and 3-programmable positions are standard
- Arms rotate for easy left or right patient entry
- Synchronized back and toe board recline
- Uphostery is designed for easy and quick clean-up
- Heavy vinyl scuff guard for toe board
- Sturdy all-metal construction
- All pivot points have permanently lubricated bearings
- Hydraulic lift system
- Quiet and dependable operation
- Chair functions are digitally controlled by a microprocessor
- Automatic return to exit
- Instant override on automatic functions
- Chair is designed for maximum serviceability
- Safety shut-off on base
- Precision photo-electric positioning

• Double articulating headrest

Articulating dual axis headrest re-positions the patient oral cavity to view the upper and lower arches, without the need of changing backrest position.



• Armrests

Offer optimum patient comfort while allowing ease of access from all directions. Armrests are movable, and can be removed completely if necessary.



• Seat Rotation Lock

Lever locks seat into desired position.



• Easy Wipe Down Contour Upholstery

Re-enforced seat upholstery provides extra comfort and high durability.

• Foot Control

Hands-free positioning the patient. 3 programmable positions plus an entry/exit position, last position and a manual disc-type positioning control.




SPECIFICATIONS

Lifting/lowering mechanism:	Hydraulic, cantilever type
Initial height/stroke:	15-1/2" (395mm) / 15" (380mm)
Backrest Reclining:	Axis mechanism
Backrest Reclining angle:	0° (horizon) – 70°
Seat tilting:	Synchronized with backrest
Seat tilting:	10° – 20°
Preset Positions:	Microcomputer controlled P1, P2, LP (last position) & AR (auto-return)
Chair rotation angle:	30° each (left and right)
Headrest:	Dual axis type
Safety switch:	Under main link
Emergency stop:	Any switch can cancel the automatic movement
Packing dimensions:	
Chair	55" x 33" x 33" / 350 lbs.
Cushion	50" x 8" x 22" / 20 lbs.

• Solid All Metal Construction

• Quiet Hydraulic System with Stable Vertical Lift

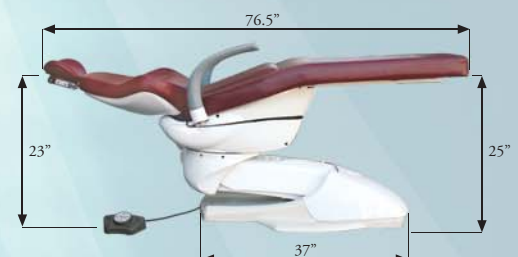

TPC

TPC Advanced Technology, Inc.
17588 E. Rowland St. • City of Industry, CA 91748

tpcdental.com

P (800) 560-8222