


— BUILT-IN CURING LIGHT SYSTEM · ALED55BN

LED 55BN

A 5 W LED curing light that installs directly in the dental unit, so it is always in the holder, always charged and always ready. Available for 110–240 V or 24 V units.

ALED55BN shown

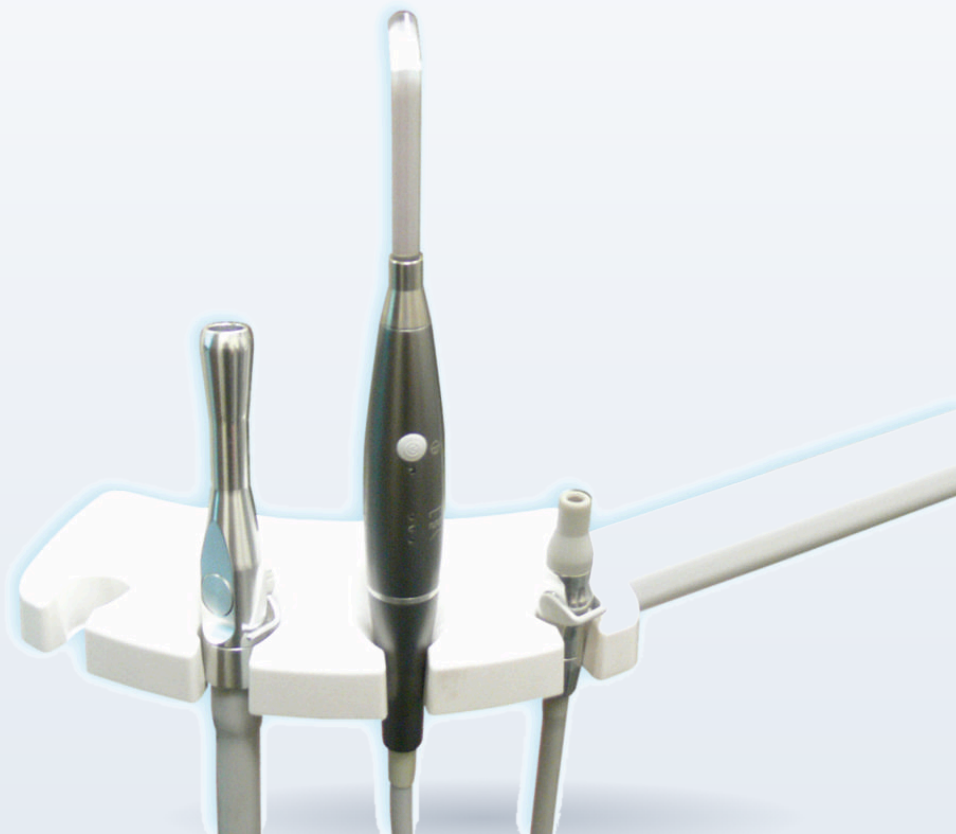

1,250 mW/cm²

Light intensity

5 W

LED light source

2 versions

110–240 V or 24 V power

2 year

Limited warranty

STANDARD FEATURES

- ✓ Installs in the dental unit
- ✓ Fast, pulse and step curing modes
- ✓ Over-temperature indicator
- ✓ No batteries or charging base
- ✓ Turbo light guide included
- ✓ Protective eye-shield included

Specifications

01 CURING MODES

Fast

F



Full power for a 20-second cycle, with a beep every 10 seconds.

Pulse

P



Interrupted output to manage heat, with a beep every 10 seconds.

Step

S



Soft start from zero to full power over 10 seconds, then full power for 10 seconds.

02 MODELS

ALED55BN-A Built-in curing light, 110–240 V

ALED55BN-E Built-in curing light, 24 V

03 SPECIFICATIONS

Light source **5 W LED**

Wavelength **430–490 nm**

Intensity **1,250 mW/cm²**

Warranty **2-year limited**

04 OPTIONS

P109-LG Light guide for LED 55B

P208 Protective eye-shield

LED 55BN-PS110 Power supply, 110–240 V

LED 55BN-PS24 Power supply, 24 V

P117 11 mm light guide

LM300 Curing light meter

05 PACKAGE INCLUDES

- ✓ Turbo light guide
- ✓ Protective eye-shield

06 PACKING DIMENSIONS

Carton	L × W × H in	Weight
Single unit	10 × 6 × 2	2 lb
Case of 20	21 × 11 × 12	27 lb

Status lights at a glance. Green means normal operation; a flashing red light means the handpiece is over temperature and needs a moment to cool.