

#### Specifications

- Battery: Lithium Ion, DC 3.6-4.2V
- Power adapter: 90-250V AC, 50-60 Hz
- Class: Class II (ICE601-1)
- Charging base: 9V DC
- Power consumption: less than 10W
- Wavelength range: 430-490 nm
- Light Intensity: 1,250 mW / cm<sup>2</sup>
- Dimensions: 7”L x 1”W x 1”H
- Weight: 160 g
- Warranty: 1 year

#### Benefits

- Ergonomic, light weight design
- High intensity output up to 1 hour before recharging
- Cordless for ultimate convenience
- Cures most materials in 5-10 seconds
- Curing time: 20 seconds per cycle with beep every 10 seconds

#### Charging Stand Light Indicators

LED color	Light Strength
Green	Normal
Orange	Light strength insufficient

#### Curing Status Light Indicators

LED color	Status Indication
Green	Normal
Red	Battery Voltage Low
Red Flashing	LED curing light handpiece is over temperature



Charging Stand

#### Curing Mode Light Indicators




LED color	Working Mode
Green	Fast curing: Time setting 20 seconds (every 10 second, it sounds once)
Red	Pulse curing: interrupt radiating (every 10 seconds, it sounds once)
Orange	Step curing: soft start for 10 second from 0 to full power, then full power during the later 10 second



Curing Light




Model LED 55N

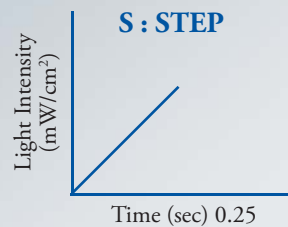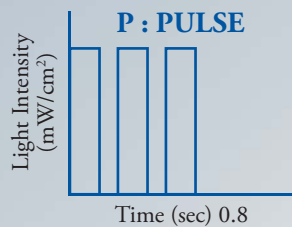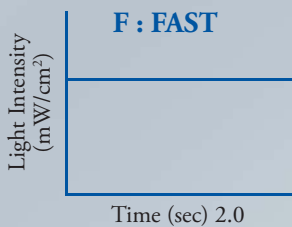

Shipping Dimensions: 10”L x 9”W x 5”H / 3 lbs.  
 (10/cs.) 25”L x 11”W x 17”H / 25 lbs.

#### Specifications

- Power adapter: 90-250V AC, 50-60 Hz
- Class: Class II (ICE601-1)
- Power consumption: less than 10W
- Wavelength range: 430-490 nm
- Light Intensity: 1,250 mW / cm<sup>2</sup>
- Dimensions: 7”L x 1”W x 1”H
- Weight: 160 g
- Warranty: 1 year

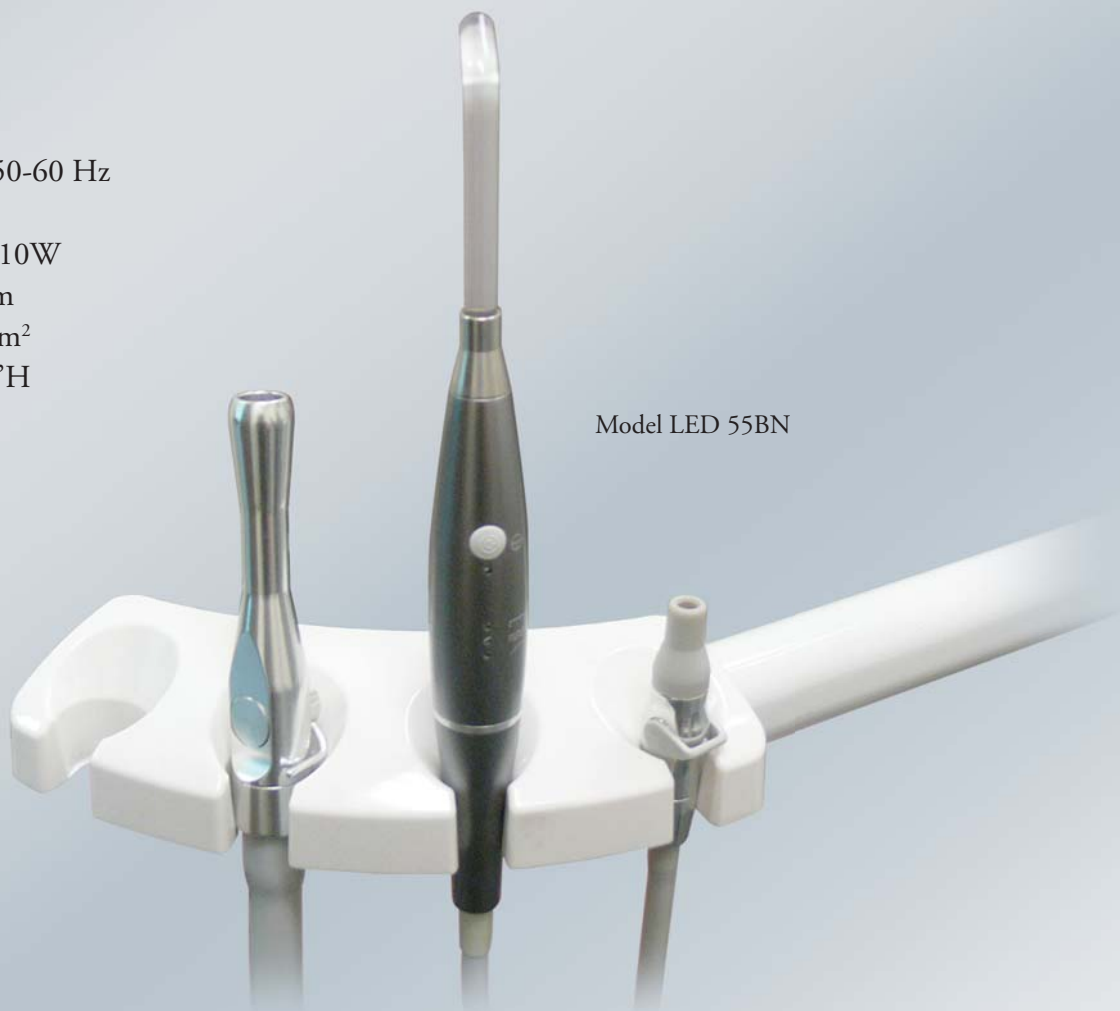




Model LED 55BN

#### Curing Mode Light Indicators

LED color	Working Mode
Green	Fast curing: Time setting 20 seconds (every 10 second, it sounds once)
Red	Pulse curing: interrupt radiating (every 10 seconds, it sounds once)
Orange	Step curing: soft start for 10 second from 0 to full power, then full power during the later 10 second

#### Curing Status Light Indicators

LED color	Status Indication
Green	Normal
Red Flashing	LED curing light handpiece is over temperature

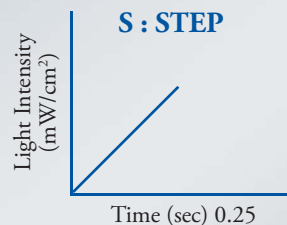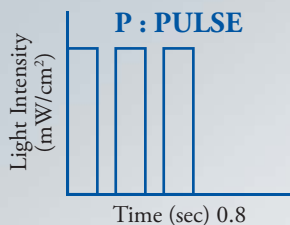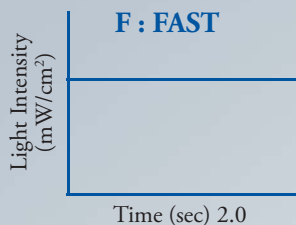

Shipping Dimensions: 10”L x 9”W x 5”H / 3 lbs.